SLENDERTONE®

Time to look and feel good

Arm Toner

INSTRUCTION MANUAL

MANUAL DE USARIO

MODE D'EMPLOI

www.slendertone.com

CONTRA-INDICATIONS, WARNINGS & PRECAUTIONS

Your SLENDERTONE product is intended to stimulate healthy muscles in order to improve or facilitate muscle performance of the triceps and biceps muscles. Your SLENDERTONE product is suitable for use by all healthy adults; however, as with other forms of exercise, some care is needed when using it. Always follow the guidelines below and read the instruction manual before use. Some of the points below are gender-specific.

CONTRA-INDICATIONS:

Do not use if:

- You have an electronic implant (e.g. cardiac pacemaker or defibrillator - as your SLENDERTONE product may interfere with the proper functioning of the implanted stimulator) or suffer from any other heart problem.
- You have any other form of medical implant. The safety has not been established for the use of electronic muscle stimulators over medical implants.
- You are pregnant. The safety of electronic muscle stimulation (EMS) over the pregnant uterus has not been established.
- You suffer from cancer. The effects of electronic stimulation on cancerous tissue are unknown.
- You are under medical supervision for cognitive dysfunction as you may not be able to comply with safety instructions.
- The unit is in close proximity (e.g., 1 yard) to shortwave or microwave therapy equipment or you are connected to high-frequency surgical equipment (e.g. diathermy), because of the risk of device malfunction due to electronic interference.
- Wearing the device necessitates placement over areas at which drugs/medicines are administered (short-term or long-term) by injection (e.g., hormone treatments).

WARNINGS:

- The long-term effects of chronic electronic stimulation are unknown.
- Unit should not be used while sleeping, because you cannot react to safety information.

When applying the GelPads™ and unit, always remember to:

- Place the GelPads and unit ONLY as indicated in this manual.

- Not place the GelPads or unit on the front or sides of the neck, across or through the heart (i.e., one GelPad on the front of the chest and one on the back), in the genital region nor on the head, because of the risk of stimulating inappropriate muscles and organs.
- Avoid any recent scars, broken or inflamed skin, areas of infection or susceptibility to acne, thrombosis or other vascular problems (e.g., varicose veins), or any parts of the body where feeling is limited.
- Avoid areas of injury or restricted movement (e.g., fractures or sprains).
- Avoid placing the GelPads over metal implants.
- Avoid application of GelPads near the thorax as it may increase the risk of cardiac fibrillation.

PRECAUTIONS:

Please read these precautions before using your SLENDERTONE unit:

- You should only use this product while in a relaxed, seated position, with your arms bent slightly and resting in your lap.
- Do not straighten your arms during the session, as this can lead to an ineffective session and may cause some discomfort.
- Never use SLENDERTONE while involved in any activity which includes using your arms or hands.
- Your arms are likely to move involuntarily during a session.
- Keep your unit in your hands during each session.
- Do not fasten the straps too tightly around your arms.
- If the garment leads are damaged or exposed, they should be replaced immediately.
- Maximum use of a SLENDERTONE product is one session per day per body area.
- Caution should be used for patients with suspected or diagnosed epilepsy.
- Caution should be used in the following:
 - When there is a tendency to bleed internally following an injury or fracture.
 - Following recent surgical procedures when electronic stimulation may disrupt the healing process
 - Over areas of skin which lack normal sensation
- Your SLENDERTONE product is not intended for medical use, for the treatment of any medical condition, nor for any permanent physical changes.
- You should only use SLENDERTONE GelPads as other kinds may reduce the effectiveness of the exercise.

Get your doctor's permission before using your SLENDERTONE product if:

- You have any serious illness or injury not mentioned in this guide.
- You have recently undergone a surgical procedure.
- You take insulin for diabetes.
- You want to use it on a young child.
- You suffer from muscle or joint problems.
- You use the unit as part of a rehabilitation program.

Adverse reactions:

- Skin irritations and burns have been reported with the use of powered muscle stimulators.
- On very rare occasions, first-time users of EMS have reported feeling light-headed or faint. We recommend that you use the product while seated until you become accustomed to the sensation.

To reposition the GelPads during a session:

- Always pause the program currently running, unfasten the strap and then refasten it around your arm after the GelPads' positions have been adjusted.

After strenuous exercise or exertion:

- Always use a lower intensity to avoid muscle fatigue.

Contact SLENDERTONE or an authorized dealer if:

- Your unit is not working correctly. Do not use in the meantime.
- You experience any irritation, skin reaction, hypersensitivity, or other adverse reaction. You should, however, note that some reddening of the skin under the garment is quite normal during and for a short time after an exercise session.

Note:

- An effective session should not cause undue discomfort.

Important:

- For first time users, muscle stimulation can be an unusual sensation. We recommend that you begin in a seated position with low stimulation intensity settings to familiarize yourself with the sensation before progressing to higher intensity settings.
- Keep your unit out of the reach of children.
- The studs and GelPads must not be connected to other objects.
- Do not use your unit at the same time as any other device which transfers an electrical current into the body (e.g., another muscle stimulator).
- Do not touch the GelPads or metal studs while the unit is switched on.
- Do not use while driving, operating machinery or cycling.
- The product is for single person use only.
- Cease using your unit if you are feeling light headed or faint. Consult your doctor if this happens
- Do not over exert yourself when using muscle stimulation. Any workout should be at a level comfortable for you.
- Although compliant with applicable EMC requirements, this device may still interfere with more sensitive equipment, please move away or turn off.
- SLENDERTONE will not accept responsibility if the guidelines and instructions supplied with this unit are not followed.

Note: If you are in any doubt about using your SLENDERTONE product for any reason, please consult your doctor before use.

WELCOME

Congratulations on purchasing your ARM TONER from SLENDERTONE and taking the first step towards firmer, stronger and more defined arms. As a doctor and the CEO of Bio-Medical Research Ltd., a company with over 40 years experience in the medical device industry, I am proud to offer you this product. You are now on the way to the body shape you've always dreamed of, so don't waste another minute - set-up your ARM TONER unit from SLENDERTONE and enjoy your first session right now!

Patricia A. Smith

Patricia A. Smith
CEO Slendertone

WHY YOU'LL LOVE IT

At SLENDERTONE we are proud to have delivered a clinically proven arm product that delivers results. This interchangeable system lets you attach the same hand-held controller to a whole range of garments, which firm, strengthen and tone different areas of your body. See www.slendertone.com for more details. The controller can be shared with your friends or partner when they purchase a SLENDERTONE SYSTEM accessory product allowing you to build and customize a complete body training collection.

YOUR ARM TONER from SLENDERTONE uses clinically proven and patented EMS Technology. A signal is sent to the nerves, through the GelPads, which in turn contracts the biceps and triceps muscles.

ARM TONER from SLENDERTONE is an electronic muscle stimulator intended to stimulate healthy muscles in order to improve or facilitate muscle performance of the triceps and biceps muscles of the upper arms.

MAKING IT WORK FOR YOU

To get the most from your ARM TONER, follow the 30 day plan and use your unit 3 times a week at an intensity level comfortable for you. We recommend you leave at least 48 hours between sessions. Use as part of a normal, healthy lifestyle.

Read all instructions fully before use!
Patricia A. Smith


THE SCIENCE BEHIND SLENDERTONE

YOUR ARM TONER USES SLENDERTONE'S patented C.S.I. Technology to deliver signals directly to your neuro-muscular system, causing contractions which exercise and tone your muscles.


Step 1: Signals are sent through the GelPads. The signals activate the nerves which control the biceps and triceps muscles.

Step 2: As the signal strength increases, the nerves, in turn, cause your muscles to contract and relax rhythmically.

The Result: A deep, effective muscle workout. It's like SLENDERTONE does the hard work for you!

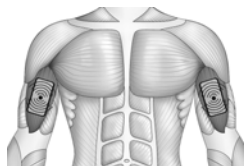BICEPS

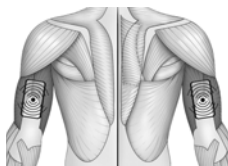
TRICEPS

WHAT'S IN THE BOX

These are the contents for the full ARM TONER pack. The Accessory Pack does not contain a unit or charger.

1. ARM TONER Unit
2. Arm Garments.
3. Medical Grade GelPads™.
4. Battery Charger
5. Carry Pouch.
6. Instruction Literature.

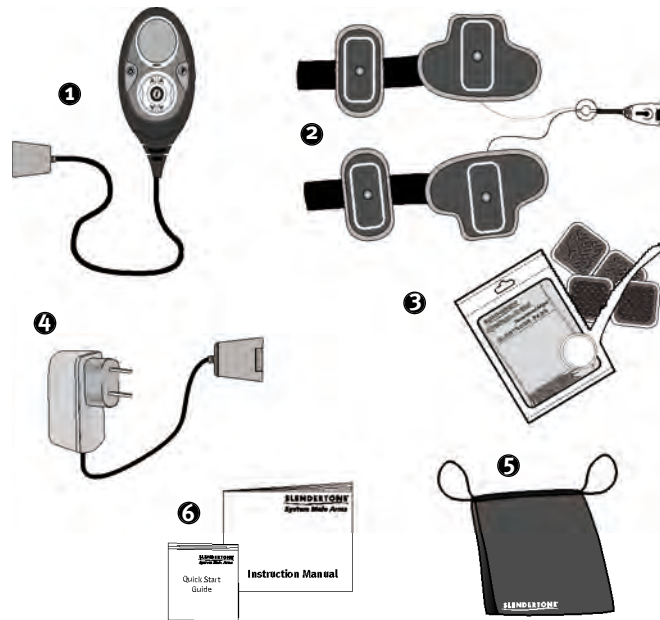UNIT CONTROLS

1. On/Off Button (⏻)

Press and hold this button for 2 seconds to turn your unit on or off. You may also pause the exercise by briefly pressing this button.

2. Increase Training Intensity (▲)

Press and hold these buttons to increase the intensity of the exercise. The right-hand button controls the garment on your right arm. The left-hand button controls the garment on your left arm.

3. Decrease Training Intensity (▼)

Press and hold these buttons to decrease the intensity of the signal to either arm.

4. Program Button (P)

Press this button to select the program you wish to use. There are 3 programs, which are explained later in this manual (page 9).

5. Information Button (i)

Press this button to see your training information, such as the maximum intensity for your current or last session and the number of sessions you have completed.





UNIT DISPLAY

- (50 Training intensity from the left arm garment.
- 46) Training intensity from the right arm garment.
- 10:08 Counts down the time left in the current session. Also displays error messages (see page 10).
- || Appears when the exercise has been paused.
- 🔇 Indicates that the sound is off.
- 🔑 Keylock is active - you cannot change the intensity levels or the program accidentally.
- ⚠ This symbol indicates a poor contact between the unit and garments or between the GelPads and the skin (see page 11 for more details).
- 🔋 Shows the battery power remaining.
- # Number of sessions completed.
- 2 Shows the program currently running (1 - 3).
- This symbol appears during each contraction phase.
- (92 87) Maximum intensity for the current/ last session.

CHARGE YOUR UNIT

1. Insert your charger into a socket and then connect it to your unit. It will take approximately 2-3 hours to fully charge the battery. Your unit is fully charged when all three sections of the battery icon are full. The battery icon will flash when the battery is low and needs to be recharged.
2. Always disconnect your unit from the charger when the battery is fully charged.

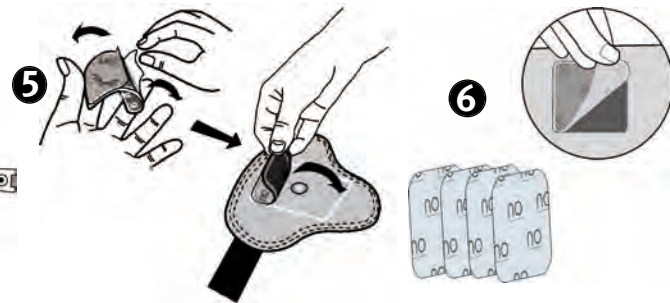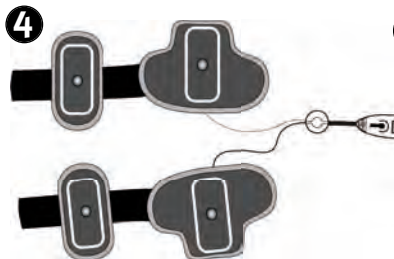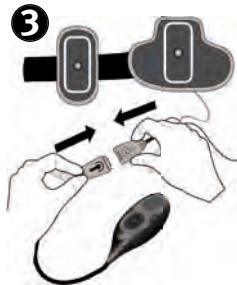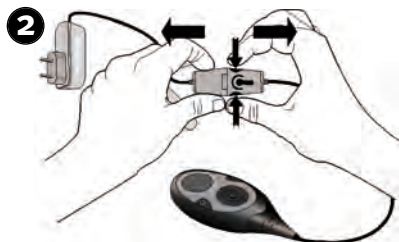**GARMENT SET-UP**


3. Connect the unit to the garments. The unit and garments must be connected properly for YOUR ARM TONER to work. "Click" the unit connector into the garment connector.
4. Place the garments with the four silver studs facing towards you. Remove the adhesive GelPads from their pack. You will notice that one side of each GelPad has a grid pattern on it while the other side is plain black.
5. Remove the covers from the **patterned** side of the four GelPads. Place the **patterned** side of the GelPads over the silver studs and within the white lines.
6. Remove the covers from the **black** side of all four GelPads. Do not throw these away as you will need to put them back onto the GelPads at the end of your session.

Note:

You should only use SLENDERTONE GelPads as other kinds may reduce the effectiveness of the exercise.

POSITIONING THE GARMENTS CORRECTLY

1. Take the garment connected to the **GOLD** lead and place the larger section centrally over the triceps muscle of your **RIGHT** arm. Ensure the wider part is towards your elbow.
2. Then, bring the strap around the **OUTSIDE** of your arm and place the small, oval section centrally over your biceps muscle. Then fasten the strap securely, but comfortably, around your arm.
3. Repeat steps 1 and 2, placing the garment connected to the **BLACK** lead onto your **LEFT** arm.

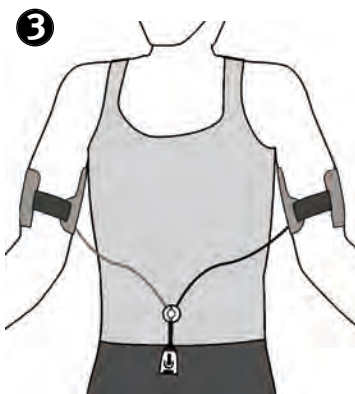4. Switch on your control unit by pressing the On/ Off button for 2 seconds. Select the program you wish to use and increase the intensity to a strong, comfortable level.

IMPORTANT:

- Before you increase the intensity, make sure that all of the silver studs are completely covered by the GelPads. Press the edges of each one firmly onto the material before use.
- You may feel a tingling sensation in your lower arms during a session. Simply changing the position of your arms can reduce this. If it persists, pause the unit and adjust the position of the arm garments slightly.

When you have completed your session, remove the garments, replace the liners on the GelPads and store the unit safely until your next session.

Important:

Always ensure the unit is switched off before removing the garments.

Have any further questions? Just visit the website: www.slendertone.com or call our Toll Free number: 1-800-551-2443

TRACK YOUR PROGRESS

Record the highest intensity levels in your performance plan (see back of this instruction manual). The diary below shows one persons plan. Each box shows the highest training intensity for each arm garment. The intensity level will vary from person to person. Fill in the diary each time you complete a training session. Leave your diary somewhere you will see it regularly. This will help to motivate you to reach your goal.

PERFORMANCE PLAN

Week 1	Session 1	Session 2	Session 3
	LEFT ARM RIGHT ARM		
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

ARM TONER PROGRAMS

YOUR ARM TONER FROM SLENDERTONE **has three programs.**

- 1. Conditioning** is a moderate strength program, introducing the biceps and triceps muscles to this form of training. This should be used before progressing to the stronger programs.
- 2. Strength** is a stronger program and will strengthen your biceps and triceps muscles as well as improving muscle performance.
- 3. Power** is a longer power program, using short, powerful contractions to increase both maximum muscle strength and sustained performance of the biceps and triceps muscles simultaneously.

- Once you have selected a program, the unit will continue using that program until you change it.
- We recommend you use the product no more than three times a week and that you leave at least 48 hours between sessions.

Note: You cannot change a program during a session, you must switch your unit off and then on again. You can then select a different program by pressing the program button.

Program	Duration	Intensity
Conditioning	10 minutes	0-99
Strength	12 minutes	0-99
Power	15 minutes	0-99

ADDITIONAL UNIT FUNCTIONS

Information Button (i)

Press the information button any time during a session to see the highest intensity level you have reached for that session (Fig. a). Pressing the information button twice shows you the total number of sessions you have completed. While pressing the information button three times identifies the product type you are using.

Fig. a

Mute Function (M)

If you want to switch off your unit's sound effects, press and hold the program button for two seconds (Fig. b). The mute function remains active indefinitely unless manually changed. Deactivate the mute function by again pressing the program button for two seconds

Fig. b

Keylock Function (K)

If you find a comfortable exercise intensity, press and hold the information button for two seconds to lock that intensity (Fig. c). This function only remains active during the session in which it is activated. If you wish to increase the intensity further, deactivate the keylock function by again pressing the information button for two seconds.

Fig. c


Error Messages

In the unlikely event of your unit developing a problem, "Err" will appear on your unit display (Fig. d). If this occurs you should switch the unit off and then switch it back on again. It should now operate properly. If the problem persists, please call your local careline for further assistance:

Tel: 1-800-551-2443

Email: info@slendertoneusa.com

Fig. d

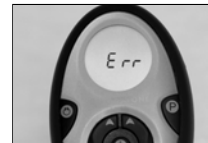Battery Power / Replacing the battery

The battery icon on the display will flash when the battery is low and needs to be recharged. After a period of time, you may find that your battery will only power your unit for a limited number of sessions, requiring you to charge the battery more frequently. This is because rechargeable batteries have a limited number of charge cycles and may eventually need to be replaced.

Should you need to replace the rechargeable battery, remove the small rubber plug from the back cover of your unit (Fig. e), unscrew the back cover and remove it. Replace the existing battery pack with a new battery pack (Fig. f) and replace the battery cover. You can purchase a new battery pack by visiting the website (www.slendertone.com) or contacting the SLENDERTONE Careline (1-800-551-2443).

Fig. e

Fig. f


TROUBLE SHOOTING CHECKLIST

If you have problems getting your ARM TONER to work, you should go through the following checklist:

Have you charged your battery?

It takes 2-3 hours to fully charge your battery

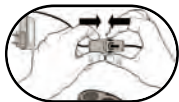Is your charger fully connected?

Push the charger connector and unit connector together until they “click” into place.

Is your unit connected to the garment properly?

Check to ensure that it is connected properly. To remove, press the outer buttons on the connector and reinsert again.

When did you last change your GelPads?

GelPads need to be changed every 20-30 sessions. Visit the website or call our support line to order replacement GelPads.

Have you removed all the plastic covers from both the grid and black sides of the GelPads?

Remember to keep one set of covers to protect the black side of the GelPads during storage.

Are the GelPads covering the metal studs?

Double-check the GelPads on your belt to ensure that the grid side of each pad is fully covering the metal studs.

Has the Pad Contact symbol (⚠) appeared on the display? *

The unit and garments are not properly connected. The ⚠ symbol always appears in conjunction with the left/ right intensity indicators. If the left indicator appears, there is a problem with the left-arm garment, the right indicator acknowledges a problem with the right-arm garment. If both indicators appear there is a problem with both garments.

Still having trouble? Visit our website or call the toll free careline:

Web: www.slendertone.com

Tel: 1-800-551-2443

Email: info@slendertoneusa.com

CARING FOR YOUR GARMENTS

Your garments may be cleaned using a lightly dampened sponge, but you must first remove the GelPads. Always follow the cleaning instructions when cleaning the garments.

Never machine-wash or hand-wash your garments. Always clean them using a lightly dampened cloth or sponge.

Do not use bleach when washing the garments.

Do not dry clean your garments.

Do not tumble dry your garments. Dry them on a flat surface. Do not dry them over anything hot. (e.g. a radiator) as they contain plastic parts. Ensure the garments are completely dry before using them again.

The garments should never be ironed.

FREQUENTLY ASKED QUESTIONS

The stimulation is uncomfortable. How can I improve this?

- Make sure the GelPads are positioned correctly and that they are pressed firmly against your skin - see pages 7-8. Switch off your unit and reposition the garments if necessary.
- You can also smear a few drops of water on the black surface of the GelPads. This can improve the comfort of the muscle stimulation, but be careful not to get water on the unit. Ensure the unit is switched OFF before you do this!
- Ensure the metal studs are fully covered by the GelPads.

My skin is red after the exercise. Is this a problem?

- Some redness of the skin after a training session is normal. It is partly due to an increase in the blood flow under the skin and should fade after a while. You may also experience some reddening of the skin due to the pressure of the garments. This is the same as the pressure marks you can get from tight clothing. You should not be concerned about this. It should fade soon after you remove them.
- If the redness is excessive, you may have the training intensity too high. This may increase the reddening in sensitive skin. Try using a lower intensity for a few days. If the problem persists, you should stop using the unit.

Will the product cause muscle soreness?

- As with all exercise some muscle soreness can occur after using your ARM TONER. This is normal and should go away after a while. Use a low intensity for a few sessions if you are experiencing muscle soreness.

How do I know when to replace the GelPads?

- With time the GelPads pick up skin debris and may need to be replaced as this makes the workout less effective and less comfortable.
- You may notice the signal is weakening even if the batteries are OK. This usually indicates that the GelPads are wearing and will soon need replacing. You can check this by first adjusting the garments to ensure correct positioning of the GelPads. If the problem persists, try using the unit with a new battery pack. If the signal is still weak or uncomfortable, you should order new GelPads.
- New GelPads can be purchased online at www.slendertone.com (check for special offers) or from the careline.

Battery performance has degraded significantly.

- After an extended period of time you may notice a degradation in the performance of the battery in your unit. At this point you should purchase a new rechargeable battery pack. New rechargeable battery packs can be purchased by visiting the website at www.slendertone.com or by calling the SLENDERTONE Careline on 1-800-551-2443.

PRODUCT WARRANTY

Dear Customer,

Thank you for purchasing your Slendertone product ("Product"). To be covered by this Limited Warranty, the Product must have been purchased directly from Slendertone Distribution, Inc. ("SDI"), or its authorized retailer in the United States. This Limited Warranty extends only to the original purchaser and cannot be transferred to anyone. Please read this carefully before using the Product.

Our Limited Warranty to you

SDI WARRANTS TO YOU, THE ORIGINAL PURCHASER, THAT UNDER NORMAL USE AND SERVICE, THE PRODUCT WILL BE FREE FROM DEFECTS IN MATERIAL AND WORKMANSHIP FOR A PERIOD OF ONE (1) YEAR FROM THE DATE OF PURCHASE ("WARRANTY PERIOD"). IF YOU DISCOVER A DEFECT COVERED BY THIS LIMITED WARRANTY DURING THE WARRANTY PERIOD, AND RETURN THE PRODUCT TO US, WE WILL REPAIR OR REPLACE THE PRODUCT, AT OUR OPTION, USING NEW OR REMANUFACTURED COMPONENTS. THERE IS NO WARRANTY OR CONDITION OF ANY KIND AS TO ANY DEFECTS DISCOVERED AFTER THE WARRANTY PERIOD EXPIRES.

Certain Product failures not covered

This Limited Warranty covers defects in manufacture that arise from correct use of the Product per the product manual and instructions. It does not cover damage caused by abuse, misuse, accidents, modification or repair by any person other than SDI, moisture, extreme heat or cold, corrosive environments, high voltage surges from external sources or an unsuitable power supply. This Limited Warranty does not cover normal wear and tear on consumable items (e.g. GelPads, garments, etc.). This Limited Warranty does not apply to any battery provided with the Product, or if the Product has an altered or defaced serial number.

How to obtain service under this Limited Warranty

To obtain service under this Limited Warranty, please contact SDI's Customer Careline at 1-800-551-2443 to obtain return instructions. You must submit proof of purchase in the form of a bona fide, dated bill of sale, receipt, or invoice (or copy thereof) evidencing that your request for service is made within the Warranty Period.

Except as may otherwise be required by legislation in your State, SDI will not be responsible for the cost of returning the Product or for any damage to or loss of the Product while it is in transit to us.

We recommend that you use a carrier that will track your shipment and ensure the Product's safe delivery. SDI cannot accept any returns sent C.O.D.

THIS LIMITED WARRANTY IS THE ONLY WARRANTY FOR THE PRODUCT, AND THERE ARE NO OTHER EXPRESS WARRANTIES, ORAL OR WRITTEN, PROVIDED BY SDI. ANY IMPLIED WARRANTIES APPLICABLE TO THE PRODUCT, INCLUDING ANY WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED TO THE DURATION OF THE WARRANTY PERIOD SET FORTH ABOVE.

Limitation of Remedies

YOUR EXCLUSIVE REMEDY, IN LIEU OF ALL INDIRECT, INCIDENTAL, CONSEQUENTIAL, EXEMPLARY, SPECIAL OR OTHER DAMAGES OF ANY KIND, INCLUDING FOR NEGLIGENCE, IS LIMITED TO THE REPAIR OR REPLACEMENT OF THE PRODUCT, IF DEFECTIVE, UNDER THE TERMS AND CONDITIONS OF THIS LIMITED WARRANTY.

SDI will bear no other expenses.

SDI reserves the right to refuse to repair or replace the Product if it is returned in an unsanitary condition. SDI reserves the right to substitute a product of equal or better quality if an identical product is not available at the time of replacement. This Limited Warranty is non-transferable.

Applicability of state law

Some states do not allow limitations on the duration of implied warranties, or the exclusion or limitation of incidental or consequential damages, and so the related limitations and exclusions stated above may not apply to you.

Governing law

This Limited Warranty is valid in the United States only and is not applicable in other countries. If you acquired the Product in the United States, the laws of the State of New York will apply to this Limited Warranty.

THIS LIMITED WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS AND YOU MAY ALSO HAVE OTHER RIGHTS, WHICH VARY FROM STATE TO STATE.

TECHNICAL SPECIFICATIONS

Caring for your unit

Your unit should not be allowed to get wet or be left in excessive sunlight. It may be cleaned regularly using a soft cloth, lightly dampened in soapy water. Do not allow the interior of your unit to become wet. Do not use detergents, alcohol, spray aerosols or strong solvents on your unit.

Access to the interior of the unit is not required for maintenance purposes.

If your unit is damaged, you should not use it but should return it to SLENDERTONE or your local distributor for replacement or repair. Repairs, service and modifications may not be carried out by anyone other than qualified service personnel authorized by SLENDERTONE.

IMPORTANT: Under no circumstance should anything other than the correct type of batteries – (rechargeable batteries 3.6V NiMH) be used with your unit. These can be purchased from the SLENDERTONE careline.

Disposing of the GelPads and batteries

Used GelPads and batteries must never be disposed of in a fire but in accordance with your state's laws governing the disposal of such items.

NOTES: The signal gradually increases to a peak intensity level at the start of the contraction phase (i.e. ramp up) and gradually decreases to nothing at the start of the relaxation phase (i.e. ramp down). When the  symbol is flashing, the stimulus is significantly reduced.

Product Type: 397

Intended use: Muscle stimulator

Classification: Internally powered equipment, Class II charger, Type BF applied parts.

Waveform: Symmetrical bi-phasic square waveform when measured into a resistive load

Accessories:

You can purchase all accessories at slendertone.com. Under no circumstances should anything other than SLENDERTONE accessories be used with your SLENDERTONE charger or SLENDERTONE unit (Type 390). Any others may not be compatible with your unit and could degrade the minimum safety levels.

SLENDERTONE SYSTEM-ABS Female Garment (E:10)

SLENDERTONE SYSTEM-ABS Male Garment (X:10)

SLENDERTONE SYSTEM-ABS GelPads:

- 1 large GelPad (Type 706 or 709)
- 2 small GelPads (Type 707 or 710)

SLENDERTONE SYSTEM-BOTTOM TONER Garment (E:30)

SLENDERTONE SYSTEM-BOTTOM TONER GelPads: • 4 large GelPads (Type 716)

SLENDERTONE SYSTEM-ARMS (E:60) Garments

SLENDERTONE SYSTEM-ARMS (E:60) GelPads: • 4 small GelPads (Type 715)

SLENDERTONE ARM TONER (X:60) Garments

SLENDERTONE ARM TONER (X:60) GelPads: • 4 small GelPads (Type 718)

SLENDERTONE Belt Extension

Battery Pack (3.6V, NiMH)

Battery Charger (EU - 2504-0303)

Battery Charger (US/Japan - 2504-0302)

Charger complies to EN 60950, Input: nominal voltage 100-240V, frequency 50-60Hz, current 125mA. Do not use any other chargers or any other power supplies.

English

Environmental Specifications:

Operating Range:	Temperature: 0 - 35°C (32 - 95°F)
	Humidity: 20 - 65 % RH
Storage Range:	Temperature: 0 - 35°C (32 - 95°F)
	Humidity: 20 - 85 % RH

Description of your unit's symbols:

There are a number of technical markings on your unit. These can be explained as follows:

The unit and garments are manufactured in China for Bio-Medical Research Ltd, Parkmore Business Park West, Galway, Ireland.

The unit requires 1 x 3.6 volt (NiMH) DC battery pack. DC is indicated by the symbol: 

The Output Frequency indicates the number of pulses per second transmitted by the unit. This is measured in hertz, which is indicated by "Hz".

Output (RMSA) means the maximum output root mean square current for each channel.

Output (RMSV) means the maximum output root mean square voltage for each channel.

This symbol means "Attention, consult the accompanying documents".

This symbol means type BF applied parts.

This symbol on your unit is to indicate conformity to the requirements of the Medical Device Directive (93/ 42/EEC). 0366 is the number of the notified body (VDE).

SN stands for 'serial number'. Inside the box, on the back of the unit is the serial number specific to this unit. The letter preceding the number indicates the year of manufacture, where "L" denotes 2006, "M" denotes 2007 etc.

The garment's batch number is represented on the garment packaging, by the number corresponding with the LOT symbol.

Rated Outputs – Voltage/Currents: X:60

Parameter	500Ω	1KΩ	1K5Ω
Output RMSV	3.8V	7V	7.9V
Output RMSA	7.7mA	7mA	5.3mA
Output Frequency	50Hz	50Hz	50-Hz
DC Component: Approx.	0 C	0 C	0 C
Positive Pulse Width:	100-125µs	100-125µs	100-125µs
Negative Pulse Width:	100-125µs	100-125µs	100-125µs
Interphase Interval:	100 µs	100 µs	100 µs

Attention:

This unit can deliver currents in excess of 2mA/cm² if used with incorrect electrodes.

At the end of the product lifecycle, do not throw this product into the normal household garbage, but bring it to a collection point for the recycling of electronic equipment.

Some product materials can be re-used if you bring them to a recycling point. By re-using some parts or raw materials from used products you make an important contribution to the protection of the environment. Please contact your local authorities if you need more information about collection points in your area.

Waste Electrical and Electronic Equipment can have potentially harmful effects on the environment. Incorrect disposal can cause harmful toxins to build up in the air, water and soil and can be harmful to human health.

Garment Materials:

Outer material: 80% Nylon, 20% polyester; Inner Material: 100% Polyester laminated to 100% polyurethane; Binding: 100% Nylon; Hook & Loop: 100% Nylon; Foam: EVA; Elasticated Straps: 100% Nylon.

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TRACK YOUR PROGRESS

Record the highest intensity levels in your performance plans (see page 9). Each box should show the highest training intensity for each arm garment. The intensity level will vary from person to person. Fill in the diary each time you complete a training session. Leave your diary somewhere you will see it regularly. This will help to motivate you to reach your goal.

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Week 1	Session 1	Session 2	Session 3
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

Week 1	Session 1	Session 2	Session 3
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

Week 1	Session 1	Session 2	Session 3
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

Week 1	Session 1	Session 2	Session 3
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

Week 1	Session 1	Session 2	Session 3
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

Week 1	Session 1	Session 2	Session 3
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

Week 1	Session 1	Session 2	Session 3
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

Week 1	Session 1	Session 2	Session 3
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

Week 1	Session 1	Session 2	Session 3
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

Week 1	Session 1	Session 2	Session 3
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

Week 1	Session 1	Session 2	Session 3
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

Week 1	Session 1	Session 2	Session 3
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

Week 1	Session 1	Session 2	Session 3
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

Week 1	Session 1	Session 2	Session 3
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

Week 1	Session 1	Session 2	Session 3
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

Week 1	Session 1	Session 2	Session 3
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

Week 1	Session 1	Session 2	Session 3
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

Week 1	Session 1	Session 2	Session 3
Week 2	Session 1	Session 2	Session 3
Week 3	Session 1	Session 2	Session 3
Week 4	Session 1	Session 2	Session 3

SLENDERTONE[®]

Arm Toner

1-800-551-2443

info@slendertoneusa.com

Designed by & Manufactured for:

Bio-Medical Research Ltd., Parkmore Business Park West, Galway, Ireland.

Copyright: © 2011, Bio-Medical Research Ltd. All Rights Reserved

Part No: 2400-0401 Rev.: 2 Date of Issue: 4/11