DENTIST INSPIRED
CUPPING ACTION

Oral-B®

PROFESSIONAL
CARE

powered
by

BRAUN

SmartSeries 4000

• PRESSURE SENSOR

English	4
Español	12
Français	19

Internet:


www.oralb.com
www.braun.com

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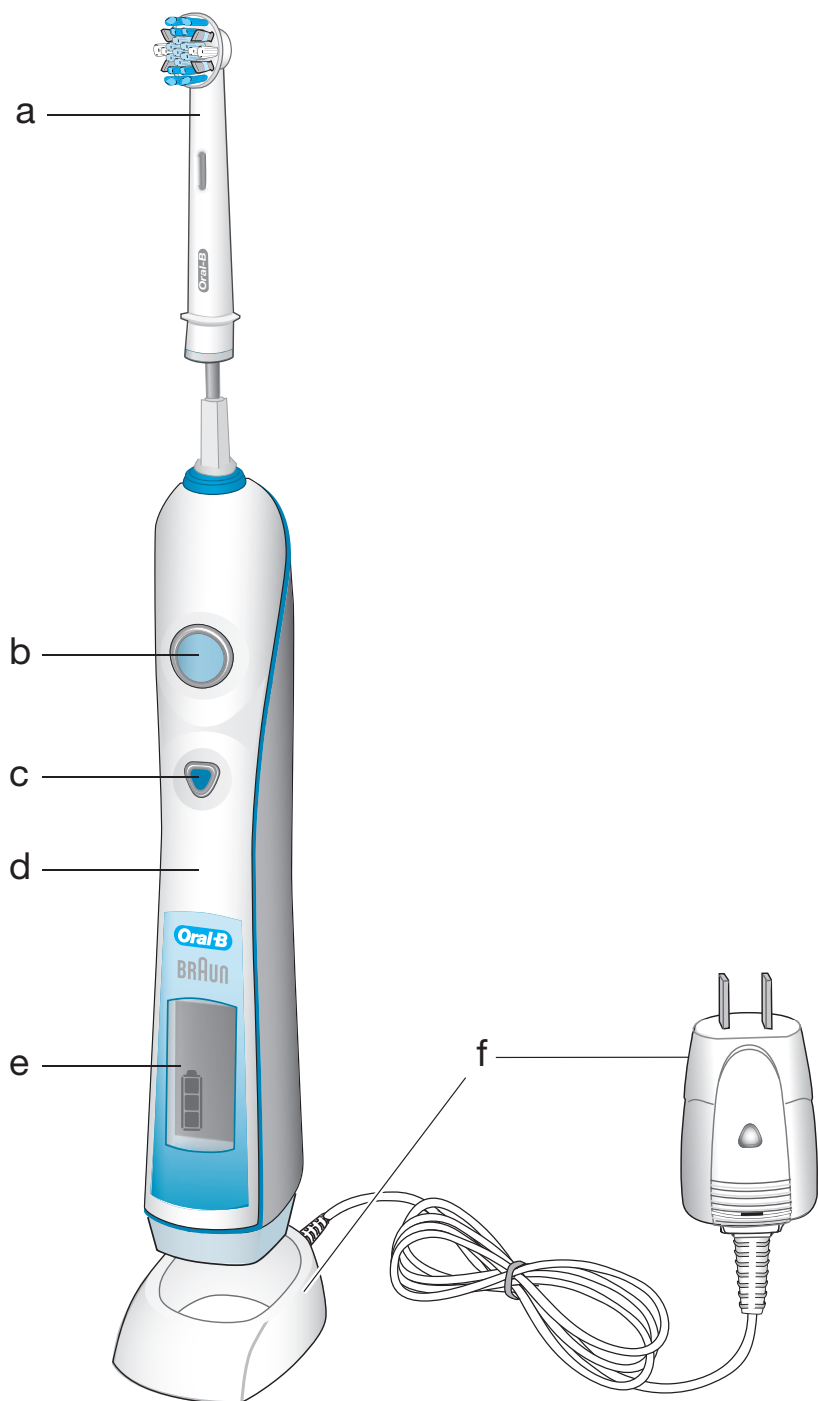IMPORTANT SAFEGUARDS

Oral-B® SmartSeries 4000® has been carefully designed to offer you and your family a unique brushing experience that is both safe and effective. When using electrical products, however, some basic safety precautions should always be followed.

Read all instructions before using.

DANGER

In order to reduce the risk of electrocution:

1. Do not place or store the charger where it can fall or be pulled into a tub or sink.
2. Do not place the charger in water or other liquid.
3. Do not reach for a charger that has fallen into water. Unplug immediately.

WARNING

To reduce the risk of burns, electrocution, fire or injury:

1. This appliance is not intended for use by children or persons with reduced physical or mental capabilities, unless they are supervised by a person responsible for their safety.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Keep cord away from heated surfaces.
4. Never drop the appliance or insert any object into any opening of the appliance.
5. Never operate this product if it has been dropped or damaged, if it is not working properly, or if the SmartPlug is damaged or has been dropped into water. If any of these occur, please contact a service center and return the product for examination and repair.
6. If the product is dropped, the brush head should be replaced before the next use even if no damage is visible.
7. This appliance is provided with a SmartPlug with integrated Safety Extra Low Voltage power supply. Do not exchange or manipulate any part of it otherwise there is the risk of electric shock. Periodically check the SmartPlug for damage. A damaged or non-functioning appliance should no longer be used. Take it to an Oral-B Braun Service Center.

IMPORTANT

1. Do not touch the power plug with wet hands. This can cause an electric shock.
2. When unplugging, always hold the power plug instead of the cord.
3. Do not modify the product. Do not disassemble or repair. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Service center. Do not disassemble the product except when disposing of the battery.
4. When taking out the battery for disposal, use caution not to short the positive (+) and negative (-) terminals.
5. This product is not intended for use by children under age three.

CAUTION

If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

SAVE THESE INSTRUCTIONS

Description

- a Brush head (depending on model)
- b On/off button
- c Brushing mode button
- d Handle
- e Handle display
- f SmartPlug charger (includes charging ring and SmartPlug)

Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.



- Plug the SmartPlug into an electrical outlet and place the toothbrush handle (d) on the charging unit.
- A battery symbol should appear in the handle display (e). The battery symbol will fill up while the handle is being charged. Once the handle is fully charged, the battery symbol will be filled with three bars. A full charge takes at least 10 hours and enables approximately 10 days of regular brushing (twice a day, 2 minutes).







Please note: If the battery symbol is not displayed immediately, continue charging. It will appear in 10–15 minutes.

- When the rechargeable battery is low, the battery symbol with one remaining bar starts flashing. Your toothbrush will slow down and then stop. It is time to recharge your toothbrush. When the battery is fully discharged, the empty frame of the battery symbol flashes. It will then need a 15-minute charge before you can use it for a 2-minute-brushing.
- For everyday use, the toothbrush handle can be stored on the charging unit to maintain it at full power. There is no risk of overcharging the battery.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

For electric specifications, see printing on the SmartPlug.

English

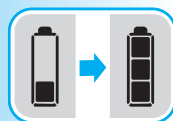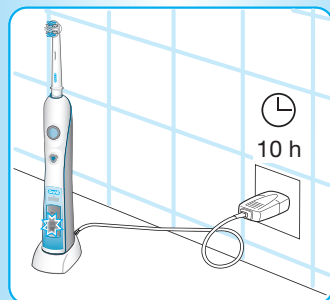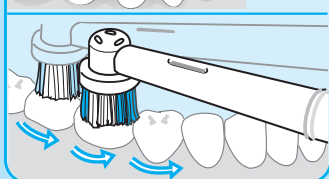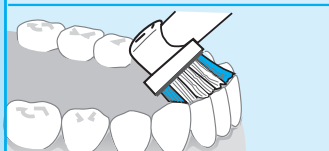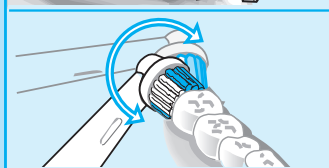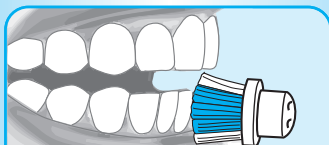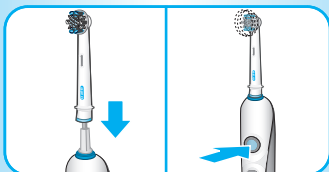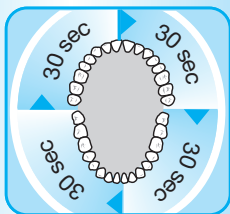Professional Timer

For an overall clean, the «Professional Timer» helps you brush the four quadrants of your mouth equally. It signals with a short stuttering sound at 30-second intervals (45-seconds intervals within «DeepClean» mode) that it is time to move to the next quadrant of your mouth, followed by a long stuttering sound when you reach the recommended 2-minute mark.

The «Professional Timer» memorizes the elapsed brushing time, even when the handle is briefly switched off during brushing. When pausing longer than 30 seconds, the timer resets to the initial starting point.

In addition, a «smiley» is displayed when you reach the dental professional recommended 2-minute brushing time. The «smiley» will be displayed for 30 seconds before it disappears automatically. Alternatively, you can directly clear the «smiley» on the display by pushing the mode button.

Brushing

You can start your toothbrush by pressing the on/off button (b). To avoid splashing, guide the brush head to your teeth before switching the appliance on.

Brushing technique

- Guide the brush head slowly from tooth to tooth.
- Hold the brush head in place for a few seconds before moving on to the next tooth.
- Brush the gums as well as the teeth, first the outsides, then the insides, finally the chewing surfaces.
- Do not press too hard or scrub, simply let the brush do all the work.

Oral-B SmartSeries 4000 has a pressure control feature to help ensure that the right degree of pressure is applied. This helps protect your gums from excessive brushing and makes your toothbrush easy to operate. If too much pressure is applied, the oscillating movement of the brush head will continue but the pulsation will stop. In addition to feeling the change, you can also hear a different sound.

In the first days of using your toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist. If you have sensitive teeth and gums, Oral-B recommends using the «Sensitive» mode.

Brushing modes (depending on model)

Your toothbrush offers different brushing modes for varying oral care needs:



- «Daily Clean» – Exceptional mouth cleaning.
- «Sensitive» – Gentle, yet thorough cleaning for sensitive areas.
- «Whitening» – Polishing and whitening for occasional or everyday use.
- «DeepClean» – Exceptional mouth cleaning with extended brushing time. Within this mode the brushing time has been extended to 45 seconds per quadrant for a total brushing time of 3 minutes.

How to switch between brushing modes:


- When pressing the on/off button (b) your toothbrush automatically starts in the «Daily Clean» mode.
- To switch to the other modes, successively press the mode button (c).
- To return from any mode to the «Daily Clean» mode, press and hold the mode button.
- To turn your brush off, press the on/off button.

During brushing, the handle display shows the brushing mode you have selected.


Please note: The «Whitening» mode is best used with an Oral-B ProWhite brush head (see «Brush heads»).

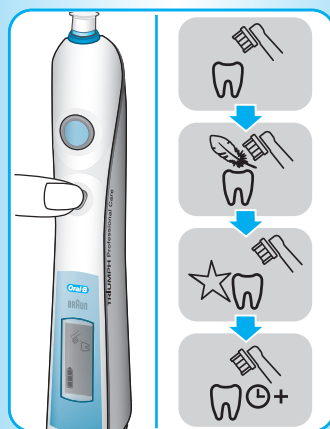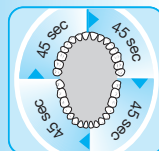Brush heads (depending on model)

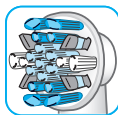The Oral-B FlossAction brush head is recommended for everyday deep cleaning. It features micropulse bristles that get deep in hard-to-reach areas between your teeth, gently removing hidden plaque.

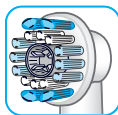For occasional or everyday whitening, please use the Oral-B ProWhite brush head. It features a dentist-inspired polishing cup to remove surface stains from coffee, tea and tobacco for effective brightening.

Please note that children under 12 years old should not use the Oral-B ProWhite brush head.



Do not use the Oral-B FlossAction or Oral-B ProWhite brush head with braces.

We recommend using the Oral-B Ortho brush head, which is specifically designed to clean around braces and wires.

The FlossAction™ and ProWhite™ brush heads feature blue Indicator® bristles helping you to monitor brush head replacement need. With proper brushing, twice a day for two minutes, the blue color will fade halfway approximately within 3 months, signalling that it is time to change the brush head.

Cleaning recommendations

- After brushing, rinse your brush head thoroughly under running water.
- Remove the brush head from the handle and clean both parts separately under running water. Then, wipe them dry before you put the handle on the charging unit.
- From time to time, you can clean the charging ring and SmartPlug with a damp cloth. The SmartPlug charger should never be placed in water or in the dishwasher.

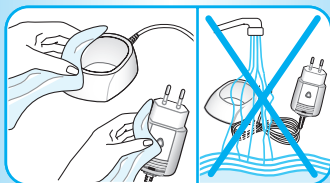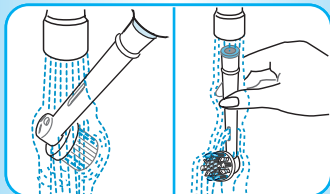Environmental notice

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Center or at approved recycling or disposing locations in accordance with local regulations.

Removing batteries

At the end of the product's useful life, open the handle as shown, remove the rechargeable battery and dispose of it according to local environmental regulations.

Caution! Opening the handle will destroy the appliance and invalidate the warranty.

Subject to change without notice.

Satisfaction Guaranteed!

Take the challenge! Experience the Oral-B difference! We guarantee that you, your dentist, and your dental hygienist will notice a difference.

Oral-B **Professional Care SmartSeries 4000** is Oral-B's most technologically advanced rechargeable toothbrush for outstanding cleaning and whitening, gentle gum care and freshening of your mouth. We are so convinced you'll fall in love with your new toothbrush that we give you this risk-free offer.

Try **Oral-B Professional Care SmartSeries 4000** for 60 days, starting from the day of purchase. If you are not 100% satisfied, receive a full refund. Guaranteed. For a full refund, please return your handle, refills and charging unit with your original cash register receipt to the address below postmarked within 60 days of purchase.


U.S. only:

Attention: 60-Day Challenge
Oral-B Braun Division
The Procter & Gamble Company
Sobin Park and A Street, MS 1F
South Boston, MA 02127

Canada only:

Attention: 60-Day Challenge
Braun Canada
a division of Procter & Gamble Inc.
Toronto, ON
M5W 1C5

For further details, please contact our customer service
1-800-566-7252.

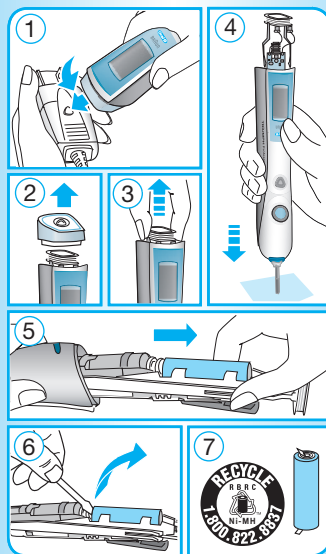FOR USA ONLY

Limited 2-Year Warranty

Oral-B Braun warrants this appliance to be free of defects in material and workmanship for a period of two years from the date of original purchase at retail.

If the appliance exhibits such a defect, Oral-B Braun will, at its option, repair or replace it without cost for parts and labor. The product must be carried in or shipped, prepaid and insured (recommended), to an authorized Oral-B Braun Service Center. Proof of date of original purchase is required.

The cost to ship the product to an authorized Service Center is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Oral-B Braun.

This warranty does not cover finishes or normal wear, nor does it cover damage resulting from accident, misuse, dirt, water, tampering, unreasonable use, service performed or attempted by unauthorized service agencies, failure to provide reasonable and necessary maintenance, including recharging, or units that have been modified or used for commercial purposes.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL ORAL-B BRAUN BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you. This warranty gives you specific legal rights, and you may have other legal rights which vary from state to state.

How to obtain service in the USA

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Braun Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Oral-B Braun Service Center of your choice or,
- B. Ship the product to the authorized Oral-B Braun Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

Declaration of conformity

This appliance complies with part 15 of the FCC rules. Operation is subject to the following two conditions: (1) this appliance may not cause harmful interference, and (2) this appliance must accept any interference received, including interference that may cause undesired operation.

Display unit: Tested to comply with FCC standards. FOR HOME OR OFFICE USE.

FOR CANADA ONLY

Limited 2-Year Warranty

In the event this appliance fails to function within the specified warranty period because of defects in material or workmanship, and the consumer returns the unit to an authorized service center, Braun Canada will, at its option either repair or replace the unit without additional charge to the consumer. This warranty does not cover any product which has been damaged by dropping, tampering, servicing performed or attempted by unauthorized service agencies, misuse or abuse. This warranty excludes units which have been modified or used for commercial purposes.

To the extent allowed by law, Braun Canada shall not be responsible for loss of the product, loss of time, inconvenience, commercial loss, special or consequential damages. The provisions of this warranty are in addition to and not a modification of or subtraction from the statutory warranties and other rights and remedies contained in applicable provincial legislation.

How to obtain service in Canada

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Braun Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Oral-B Braun Service Center of your choice or,
- B. Ship the product to the authorized Oral-B Braun Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

Declaration of conformity

Operation is subject to the following two conditions: (1) this appliance may not cause harmful interference, and (2) this appliance must accept any interference received, including interference that may cause undesired operation.

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