SALTER

MiBody

BLOOD PRESSURE MONITOR

Instructions and Guarantee

BPW-9154

Downloaded from www.Manualslib.com manuals search engine

GB

MiBODY

Salter MiBody is a range of connected devices to make weight and health management easier. Measurements are sent wireless to your mobile device and can be viewed using the free to download MiBody App. Trends can be tracked overtime and the measurements explained with healthy/recommended ranges. The results are stored on your personal mobile device only and are private unless you choose to share your achievements.
MiBody should be used in addition with and not as a substitute to professional medical guidance. If you are concerned about any aspect of your health please seek professional medical guidance.

INDICATIONS FOR USE

The Salter MiBody Blood Pressure Monitor is intended for use in measuring blood pressure and heartbeat rate with wrist circumference ranging from 13.5 cm to 21.5 cm (about 5-8.5 inches).
It is intended for adult indoor use only.

CONTRAINDICATIONS

This device is contraindicated for any female subject who is pregnant and any children under the age of 18 years old.
This device is not intended to be a diagnostic device. Contact your physician if hypertensive values are indicated.

MEASUREMENT PRINCIPLE

This product uses the Oscillometric Measuring method to detect blood pressure. Before every measurement, the unit establishes a “zero pressure” equivalent to the air pressure. Then it starts inflating the cuff, meanwhile, the unit detects pressure oscillations generated by beat-to-beat pulsatile, which is used to determine the systolic and diastolic pressure, and also pulse rate. The device also compares the longest and the shortest time intervals of detected pulse waves to mean time interval then calculates standard deviation. The device will display a warning signal with the reading to indicate the detection of irregular heartbeat when the difference of the time intervals over.

WEEE EXPLANATION

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

BATTERY DIRECTIVE

This symbol indicates that batteries must not be disposed of in the domestic waste as they contain substances which can be damaging to the environment and health. Please dispose of batteries in designated collection points.

GUARANTEE

This product is intended for domestic use only. Salter will repair or replace the product, or any part of this product, free of charge if within 2 years of the date of purchase, it can be shown to have failed through defective workmanship or materials. This guarantee covers working parts that affect the function of the scale. It does not cover cosmetic deterioration caused by fair wear and tear or damage caused by accident or misuse. Opening or taking apart the scale or its components will void the guarantee. Claims under guarantee must be supported by proof of purchase and be returned carriage paid to Salter (or local Salter appointed agent if outside the UK). Care should be taken in packing the scale so that it is not damaged while in transit. This undertaking is in addition to a consumer’s statutory rights and does not affect those rights in any way. For UK Sales and Service contact HoMedics Group Ltd, HoMedics House, Somerhill Business Park, Five Oak Green Road, Tonbridge, Kent TN11 0GP, UK. Helpline Tel No: (01732) 360783. e-mail: support@homedics.co.uk. For Ireland, please contact Petra Brand Masters, Unit J4 Maynooth Business Campus, Maynooth, Co. Kildare, Ireland. Tel +00 353 (0) 1 6510660. e-mail sales@petrabrandmasters.ie.
www.salterhousewares.com/servicecentres

Downloaded from www.Manualslib.com manuals search engine

SAFETY INFORMATION

The below signs might be in the user manual, labeling or other components. They are the requirement of standard and using.


	Symbol for "THE OPERATION GUIDE MUST BE READ"	Bluetooth®	The Bluetooth Combination Mark
CE 0120	Symbol for "COMPLIES WITH MDD 93/42/EEC REQUIREMENTS"		Symbol for "TYPE BF APPLIED PARTS"
	Symbol for "MANUFACTURER"		Symbol for "ENVIRONMENT PROTECTION - Wast electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for recycling advice"
SN	Symbol for "SERIAL NUMBER"		Symbol for "DIRECT CURRENT"

CAUTION

- Please do read this user manual carefully and thoroughly before use.
- This device is intended for adult use only.
- This device is intended for non-invasive measuring and monitoring of asterial blood pressure. It is not intended for use on extremities other than the wrist or for functions other than obtaining a blood pressure measurement.
- Do not confuse self-monitoring with self-diagnosis. This unit allows you to monitor your blood pressure. Please start or end medical treatment basing solely on physician's treatment advice.
- If you are taking medication, consult your physician to determine the most appropriate time for your measurement. Never change a prescribed medication without your physician's consent.
- This unit is not suitable for continuous monitoring during medical emergencies or operations.
- If the pressure of the cuff exceeds 40 kPa (300 mmHg), the unit will automatically deflate. Should the cuff not deflate when its pressure exceeds 40 kPa (300 mmHg), detach the cuff from the wrist and press the START/STOP button to stop inflation.
- Do not use the monitor under the conditions of strong electromagnetic field (e.g. medical RF equipment) that radiates interference signal or electrical fast transient / burst signal.
- The maximum temperature that the applied part can be achieved is 42.5°C while the environmental temperature is 40°C.
- The device is not AP/APG equipment. It is not suitable for use in the presence of a flammable anaesthetic mixture with air (or oxygen, nitrous oxide).
- Please keep the unit out of reach of infants or children, since inhalation or swallowing of small parts is dangerous or even fatal.
- Please use ACCESSORIES and detachable parts specified / authorised by MANUFACTURER. Otherwise, it may cause damage to the unit or danger to the user / patient.

Downloaded from www.Manualslib.com manuals search engine

DISPLAY GUIDE

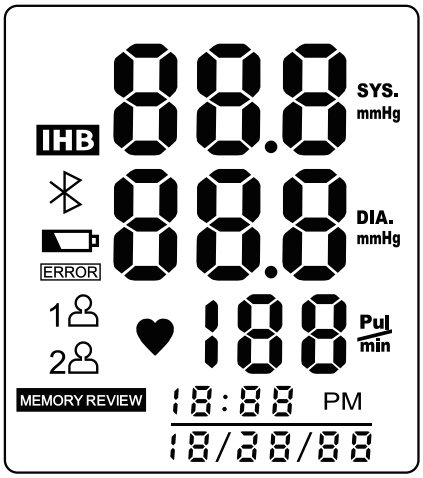SYMBOL	DESCRIPTION	EXPLANATION
SYS	Systolic Blood Pressure	
DIA	Diastolic Blood Pressure	
Pul / min	Pulse	beat/minute
+Lo	Low Battery	Low battery and please replace the batteries
mmHg	Unit	Measurement unit of blood pressure
PM	Time	Hour:Minute (Month/Day/Year)
IHB	IHB Detector	Irregular Heartbeat Detector
	Bluetooth	Successful Bluetooth Connection
ERROR	Error	Error
1 2	User ID	The selected User ID
MEMORY REVIEW	Memory	Recalling the history records

Downloaded from www.Manualslib.com manuals search engine

INSTALLING AND REPLACING THE BATTERIES

1. Open the battery door.
2. Insert the batteries according to the polarity indications.
 (Always select the authorized / specified battery: Two AAA-size alkaline batteries).
3. Close the battery door.

Replace the batteries under following circumstances:

- +Lo displays on the LCD
- The LCD display dims
- When powering on the monitor, the LCD doesn't light up.

CAUTION

- Remove batteries if the device is not likely to be used for some time.
- Worn batteries are harmful to the environment. Do not dispose with daily garbage.
- Remove the old battery from the device following your local recycling guidelines.

SETTING DATE & TIME

1. When the monitor is **off**, press and hold SET button for 3 seconds.
2. The hour digits will flash. Press MEM button to change the hour.
3. Press SET to confirm.
4. The minute digits will flash. Press MEM button to change the minute.
5. Press SET to confirm.
6. The month digits will flash. Press MEM button to change the month.
7. Press SET to confirm.
8. The day digits will flash. Press MEM button to change the day.
9. Press SET to confirm.
10. The year digits will flash. Press MEM button to change the year.
11. Press SET to confirm.
12. **dOnE** will be displayed then the monitor will turn off.

Note: To update or change the time or date repeat the above procedure making changes as required.

MiBODY APP

Before using your device for the first time:

1. Download and install the Salter MiBody App from the App store. Use keyword search terms "Salter" or "MiBody"
2. Turn Bluetooth "On" in the "Settings" menu of your iPhone/iPad.
3. Open the MiBody App and follow the on screen instructions to set up your user account.

Downloaded from www.Manualslib.com manuals search engine

PAIRING DEVICES

1. Turn Bluetooth "On" in the "Settings" menu of your iPhone/iPad.
2. Open the MiBody App.
3. Press the ⚙ button to enter "Settings".
4. Press the + button.
5. Press the ᛒ button.
6. When the blood pressure monitor is off press and hold the START/STOP button until 0 is displayed on the screen.
7. On your iPhone/iPad, touch the blood pressure monitor icon to connect with the blood pressure monitor. Select the user number you wish to pair. Touch the Done button to finish.
8. If pairing is successful [] will be displayed.
9. If pairing is unsuccessful E 1 will be displayed.
10. The blood pressure monitor will automatically switch off.
11. Repeat procedure for each user and each iPhone/iPad.

Note: Each user number can only be paired with one iOS device. Results will be sent to the paired device only.

POSITIONING THE CUFF

1. Remove all accessories (watch, bracelet, etc) from your left wrist. If your physician has diagnosed you with poor circulation in your left wrist, use your right wrist.
2. Roll or push up your sleeve to expose the skin.
3. Apply the cuff to your left wrist with your palm facing up.
4. Position the top edge of the cuff about 1-1.5cm from the wrist.
5. Fasten the wrist cuff around your wrist, leaving no extra room between the cuff and your skin. If the cuff is too loose, the measurement will not be accurate.
6. During measurement the monitor should be level with your heart.

BLOOD PRESSURE MEASUREMENT

1. After correctly positioning the cuff ensuring the monitor is level with your heart press the START/STOP button. The monitor turns on and the user number is displayed.
2. Press ▲ MEM or SET ▼ button to change the user number.
3. Press START/STOP button to confirm and start the measurement.
4. Measurement is automatic and takes approximately 30 seconds. The monitor will inflate and squeeze your wrist. After measurement the monitor will deflate and the results will be displayed on the screen.
5. If the user number is paired with an iPhone/iPad ᛒ will be displayed and the data will be sent wirelessly to your mobile device.
6. If the device is out of range or the App is closed ERROR will be displayed and the data will be saved on the monitor and transferred next time the monitor is used.

Note: To receive data the MiBody App must be open or running in the background. The MiBody App will not receive data if it is closed completely or Bluetooth is turned off.

Caution: Data is not saved or transferred for **GuESt** users.

Downloaded from www.Manualslib.com manuals search engine

RECALLING THE RECORDS

1. Press MEM button to access the memory. The monitor will display the last measurement data.
2. Press MEM button or SET button to rotate the history records.
 MEM to go forward;
 SET to go backward.

CAUTION

The most recent record (1) is shown first. Each new measurement is assigned to the first (1) record. All other records are pushed back one digit (e.g., 2 becomes 3, and so on), and the last record (60) is dropped from the list.
If you want to check the other users' measurement data. Please press the START button to turn on the monitor to select the correct user, press the STOP button to turn off the monitor. Then press "MEM" button to access the memory.
The User Guest only can check the last record.

DELETING THE RECORDS

When you did not obtain the accurate measurement, you can clear all the measuring results by following below steps.

1. Under Memory Recalling Mode, press and hold both the MEM button for 3 seconds.
2. When the LCD display blinking **dEL ALL**, press SET button to confirm clearing.
 The LCD will display **dEL dOnE**, indicating that memory clearing is complete. And then it will turn off automatically.
3. If you wish to cancel clearing, press START/STOP to turn off the monitor.
4. When there is no memory in the monitor, if you press the MEM button to look up history, the LCD will display - - -.

TIPS FOR MEASUREMENT

It can cause inaccuracy if the measurement is taken in the following circumstances.

1. Within 1 hour after eating a meal or drinking.
2. Immediate measurement after tea, coffee or smoking.
3. Within 20 minutes of taking a bath.
4. When talking or moving your fingers.
5. In a very cold environment.
6. When you need the lavatory.
7. Rest for 5 minutes before measuring
8. Wait at least 3 minutes between measurements. This allows your blood circulation to recover.
9. Try to measure under similar conditions for example, take daily measurements at approximately the same time, on the same wrist, or as directed by a physician.

Downloaded from www.Manualslib.com manuals search engine

GB

ADVICE FOR USE AND CARE

To obtain the best performance, please follow below instructions.

- Put in a dry place and avoid the sunshine.
- Avoid immersing it in the water. Clean it with a dry cloth in case.
- Avoid shaking and collision.
- Avoid dusty environment and unstable temperature surrounding.
- Use the slightly damp cloth to remove the dirt.
- Avoid washing the cuff.

Cleaning: Dust environment may affect the performance of the unit. Please use the soft cloth to remove the dirt before and after use.

Please make sure the unit functions safely and it is in proper working conditions before use.


Please follow the instructions for correct replacement of interchangeable or detachable parts specified by SERVICE PERSONNEL of MANUFACTURER as "replaceable".

Disposal: Degraded sensors may result in inaccurate measurement while loosened electrodes may cause the monitor's failure to power on. Please dispose of ACCESSORIES, detachable parts, and ME EQUIPMENT according to local guidelines.

ABOUT BLOOD PRESSURE

What are systolic pressure and diastolic pressure?

When ventricles contract and pump blood out of the heart, the blood pressure reaches its maximum value in the cycle, which is called systolic pressure. When the ventricles relax, the blood pressure reaches its minimum value in the cycle, which is called diastolic pressure.

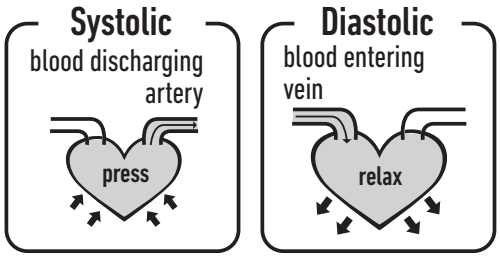
What is the standard blood pressure classification?

The blood pressure classification published by World Health Organization (WHO) and International Society of Hypertension (ISH) in 1999 is as follows:

CAUTION

Only a physician can tell your normal BP range. Please contact a physician if your measuring result falls out of the range. Kindly note that only a physician could tell whether your blood pressure value has reached a dangerous point.

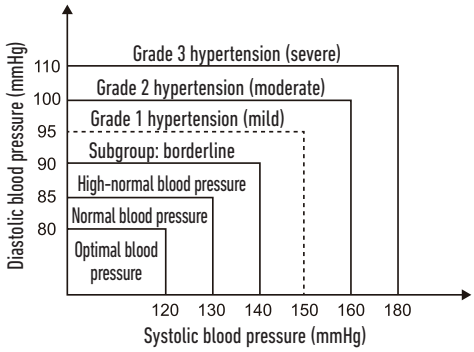Blood pressure (mm Hg) \ Level	Optimal	Normal	High-normal	Mild	Moderate	Severe
SYS	<120	120–129	130–139	140–149	160–179	≥180
DIA	<80	80–84	85–89	90–99	100–109	≥110

Downloaded from www.Manualslib.com manuals search engine

IRREGULAR HEARTBEAT DETECTOR

This Blood Pressure Monitor is equipped with an intelligent function of Irregular Heartbeat (IHB) Detector. During each measurement, this equipment records the heartbeat intervals and works out the standard deviation. If the calculated value is larger than or equal to 15, this equipment will light up the IHB symbol on the screen when displaying the measuring result.

CAUTION

The appearance of the IHB icon indicates that a pulse irregularity consistent with an irregular heartbeat was detected during measurement. Usually this is NOT a cause for concern. However, if the symbol appears often, we recommend you seek medical advice. Please note that the device does not replace a cardiac examination, but serves to detect pulse irregularities at an early stage.

QUESTIONS & ANSWERS

For a full list of MiBody frequently asked questions please visit www.uk.salterhousewares.com/mibody

Why my blood pressure varies even in one day?

1. Individual blood pressure varies during the day. It also affected by the way you tie your cuff and your measurement position, so please take the measurement under similar conditions.
2. The varies of the pressure is greater if the person take medicine.
3. Waiting at least 4-5 minutes for another measurement.

Why the blood pressure I get from the hospital is different from home?

The blood pressure is different even during 24 hour because of the weather,emotion, exercise etc, specially the “white coat” in hospital which makes the results are higher than the ones at home. The attention you need to pay when you measure your blood pressure at home:

- If the cuff is tied properly.
- If the cuff is too tight or too loose.
- If the cuff is tied on the wrist.
- If you feel anxious pressured.

You had better take deep breath 2-3 times before beginning.
Advice: adjust yourself for 4-5 minutes until you calm down.

If the result is the same if measuring on the right wrist?

It is ok for both wrists, but there will be some different results for different person, so suggest you measure the same wrist every time.

Downloaded from www.Manualslib.com manuals search engine

GB

TROUBLESHOOTING

This section includes a list of error messages and frequently asked questions for problems you may encounter with your blood pressure monitor. If the products not operating as you think it should, check here before arranging for servicing.

PROBLEM	SYMPTOM	CHECK THIS	REMEDY
No power	Display is dim or will not light up.	Batteries are exhausted.	Replace with new batteries
		Batteries are inserted incorrectly.	Insert the batteries correctly
Low batteries	+Lo Show on the display	Batteries are low.	Replace with new batteries
Error message	E 1 shows	Data communication has failed	Make sure that phone's Bluetooth is on or within the distance range
	E 2 shows	The cuff is very tight	Refasten the cuff and then measure again.
	E 3 shows	The pressure of the cuff is excess.	Relax for a moment and then measure again.
	E 9 shows	Product has not been activated.	Reactivated
	E 10 or E 11 shows	The monitor detected motion while measuring.	Movement can affect the measurement. Relax for a moment and then measure again.
	E 20 shows	The measurement process does not detect the pulse signal.	Loosen the clothing on the arm and then measure again.
	E 21 shows	Measure incorrectly.	Relax for a moment and then measure again.
	EExx, shows on the display.	A calibration error occurred.	Retake the measurement. If the problem persists, contact the retailer or our customer service department for further assistance. Refer to the warranty for contact information and return instructions.

Downloaded from www.Manualslib.com manuals search engine

SPECIFICATIONS

Power supply	2*AAA alkaline batteries
Display moder	Digital LCD V.A. 36x41mm
Measurement mode	Oscillographic testing mode
Measurement range	Pressure: 0kpa-40kpa (0mmHg-300mmHg) pulse value: (40-199) beat/minute
Accuracy	Pressure: 5°C-40°C within ± 0.4kpa (3mmHg) 0°C-45°C (out of 5°C-40°C) within ± 0.7kpa (5mmHg) pulse value: ± 5%
Normal working condition	Temperature: 5°C to 40°C Relative humidity ≤85% Atmospheric pressure: 86kPa to 106kPa
Storage & transportation condition	Temperature: -20°C to 60°C RH: 10% to 93% Atmospheric pressure: 50kPa to 106kPa
Measurement perimeter of the wrist	About 13.5cm-21.5cm
Net Weight	Approx. 120g (Excluding the dry cells)
External dimensions	Approx. 80x65x22mm (Excluding the cuff)
Attachment	2*AAA alkaline batteries, user manual
Mode of operation	Continuous operation
Degree of protection	Type BF applied part
Protection against ingress of water	IP20
Software version	V01
Device classification	Internally Powered ME Equipment

WARNING: No modification of this equipment is allowed.

Downloaded from www.Manualslib.com manuals search engine

GB

CONTACT INFORMATION

For more information about our products, please visit **www.salterhousewares.com**
For UK Sales and Service contact HoMedics Group Ltd, HoMedics House, Somerhill Business Park, Five Oak Green Road, Tonbridge, Kent TN11 0GP, UK. Helpline Tel No: (01732) 360783. For Ireland, please contact Petra Brand Masters, Unit J4 Maynooth Business Campus, Maynooth, Co. Kildare, Ireland. Tel +00 353 (0) 1 6510660. e-mail sales@petrabrandmasters.ie.
www.salterhousewares.com/servicecentres

EMC GUIDANCE

Table 1 Guidance and manufacturer's declaration – electromagnetic emissions - for all EQUIPMENT and SYSTEMS

Guidance and manufacture's declaration – electromagnetic emission		
The BPW-9154 is intended for use in the electromagnetic environment specified below. The customer or the user of the BPW-9154 should assure that it is used in such an environment.		
Emission test	**Compliance**	**Electromagnetic environment - guidance**
RF emissions CISPR 11	Group 1	The BPW-9154 must emit electro-magnetic energy in order to perform its intended function. Nearby electronic equipment may be affected.
RF emission CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Not applicable	

Downloaded from www.Manualslib.com manuals search engine

Table 2 Guidance and manufacturer's declaration – electromagnetic immunity – for all ME EQUIPMENT and ME SYSTEMS

Guidance and manufacture's declaration – electromagnetic emission			
The BPW-9154 is intended for use in the electromagnetic environment specified below. The customer or the user of the BPW-9154 should assure that it is used in such an environment.			
Immunity test	**IEC 60601 test level**	**Compliance level**	**Electromagnetic environment - guidance**
Electrostatic discharge (ESD) IEC 61000-4-2	±6 kV contact ±8 kV air	±6 kV contact ±8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines	N/A	
Surge IEC 61000-4-5	±1 kV line(s) to line(s) ±2 kV line(s) to earth	N/A	
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5% UT (>95% dip in UT) for 0.5 cycle	N/A	
	40% UT (60% dip in UT) for 5 cycles	N/A	
	70% UT (30% dip in UT) for 25 cycles	N/A	
	<5% UT (>95% dip in UT) for 5 sec	N/A	
Power frequency 3A/m (50Hz) magnetic field IEC 61000-4-8	3A/m	3A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
NOTE UT is the a.c. mains voltage prior to application of the test level.			

Downloaded from www.Manualslib.com manuals search engine

GB

Table 4 Guidance and manufacturer's declaration – electromagnetic immunity – for ME EQUIPMENT and ME SYSTEMS that are not LIFE-SUPPORTING

Guidance and manufacture's declaration – electromagnetic emission			
The BPW-9154 is intended for use in the electromagnetic environment specified below. The customer or user of the BPW-9154 should assure that it's used in such an environment.			
Immunity test	**IEC 60601 test level**	**Compliance level**	**Electromagnetic environment – guidance**
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	N/A	Portable and mobile RF communications equipment should be used no closer to any part of the BPW-9154, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. **Recommended separation distance** $d = 1.167 \sqrt{P}$
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	3 V/m	$d = 1.167 \sqrt{P}$ 80 MHz to 800 MHz $d = 2.333 \sqrt{P}$ 800 MHz to 2.5 GHz where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacture and d is the recommended separation distance in metres (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,[a] should be less than the compliance level in each frequency range.[b] Interference may occur in the vicinity of equipment marked with the following symbol:
NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies. NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			
[a] Field strengths from fixed transmitters, such as base stations for radio (cellular / cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the BPW-9154 is used exceeds the applicable RF compliance level above, the BPW-9154 should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the BPW-9154. [b] Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3V/m.			

Downloaded from www.Manualslib.com manuals search engine

Table 6 Recommended separation distances between portable and mobile RF communications equipment and the EQUIPMENT or SYSTEM – for ME EQUIPMENT or ME SYSTEM that are not LIFE-SUPPORTING

Recommended separation distances between portable and mobile RF communications equipment at the BPW-9154.			
The BPW-9154 is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the BPW-9154 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmittters) and the BPW-9154 as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter (W)	Separation distance according to frequency of transmitter (m)		
	150 kHz to 80 MHz $d = 1.167\sqrt{P}$	80 MHz to 800 MHz $d = 1.167\sqrt{P}$	800 MHz to 2.5 GHz $d = 2.333\sqrt{P}$
0.01	N/A	0.117	0.233
0.1	N/A	0.369	0.738
1	N/A	1.167	2.333
10	N/A	3.690	7.377
100	N/A	11.67	23.33
For transmitters rated at a maximum output power not listed above, the recommended separation distance *d* in metres (m) can be determined using the equation applicable to the frequency of the transmitter, where *P* is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer. NOTE 1 At 80MHz and 800MHz, the separation distance for the higher frequency range applies. NOTE 4 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

Downloaded from www.Manualslib.com manuals search engine

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