pearl by beurer

OT 80

ENGLISH Ovulation Predictor And Fertility Tracking Kit
Instructions for use 3

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Questions or Comments? Call our US-based customer service toll free at 1-800-536-0366.

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READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT

Keep this manual in a safe location for future reference

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Contents

1. Important safety notes4
2. Getting to know your instrument5
3. Package Contents5
4. Parts and Controls....6
5. To Set Up6
6. How It Works6
7. Symbols Used....7
8. To Use8
9. Care, Maintenance and Disposal....12
10. Technical Specifications13

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact Beurer customer service. Please see warranty for service contact.

Dear Valued Customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for application in the areas of heating, weight management, blood pressure, body temperature, pulse, gentle therapy, massage, and beauty. Register your product today to receive further benefits:

- added protection for your purchase
- health and wellness tips
- new product information
- continual outstanding customer service and support.

Visit www.registerbeurer.com or contact customer service at 1-800-536-0366.

Please read these instructions for use carefully and keep them for later use.

Make them accessible to other users and observe the information they contain.

With kind regards,
Your Beurer Team

1. Important safety notes

Signs and symbols

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

READ THIS ENTIRE MANUAL, THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

	This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death.
WARNING	Indicates a potentially hazardous situation, which, if not avoided, could result in death or serious injury.
CAUTION	Indicates a potentially hazardous situation, which, if not avoided, may result in minor/ moderate injury.

NOTICE	Addresses practices not related to personal injury, such as product and/or property damage.
	**Note** Note on important information
	Observe the instructions for use
	Not for re-use/for single use only
	Keep away from sunlight
	Temperature limit 39.2°F (4°C) to 86°F (30°C)
	Use by
	Batch designation
	Serial number
 / Art.-Nr.	Reference number

⚠ CAUTION - To reduce the risk of personal injury or product/property damage:

- This kit is not intended for contraceptive use.
- This kit is intended solely to predict the most fertile days of a woman's menstrual cycle. Do not use this kit for any other purpose.
- While this mobile app notifies you of the best days to conceive, that does not mean those days are the only days in which pregnancy is possible. Frequent sexual intercourse during fertile days is recommended to maximize your chances of getting pregnant.
- This mobile app is not for use when performing hormonal contraception, as hormonal contraception affects the natural course of events to achieve pregnancy.
- Use the mobile app only with the test strips included with this kit. Other test strips will produce inaccurate results.
- Using this mobile app with the included test strips is limited to countries where this kit is sold. Do not use mobile app with test strips in other countries.
- If you experience irregular cycles and use this mobile app and test strips, make sure you have enough test strips to cover the maximum length of your cycle. The kit is designed to cover cycles greater than 21 and up to 35 days. Longer cycles may require an additional set. If you have conditions, especially those that affect your natural hormone balance or cycle, please consult your doctor before using the kit.


2. Getting to know your instrument

This kit is a self-test system consisting of a mobile application that is to be used with the included FSH and LH hormonal test strips. These are designed to estimate the most fertile days of a woman's menstrual cycle. The App tracks factors of this cycle based on user input. The hormonal test strips have their own instructions for use and must be read carefully before use. These instructions for use are for this system (app and strips) only and may refer to the test strip instructions for guidance only. The App is NOT to be used for contraceptive purposes; nor alone, nor in combination with the test strips. If you are seeking natural family planning, consult a medical professional.

3. Package Contents

1 x reusable stripholder OT 80
1 x activation QR code for the Pearl Fertility app
15 x Pearl Fertility test strips - FSH
15 x Pearl Fertility test strips - LH
2 x Pearl Fertility test strips - hCG (pregnancy tests)
1 x Instructions for use (system)
1 x Instructions for use (test strips)

4. Parts and Controls

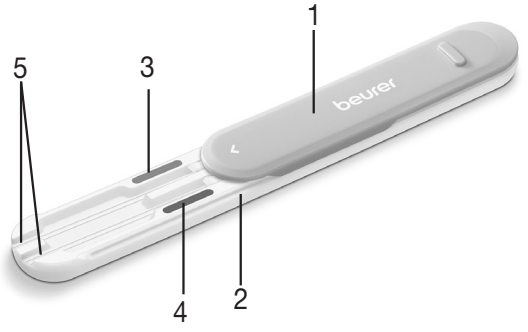
1. Removable Protective Cover
2. Test Strip Holder
3. Pink Color Indicator
4. Blue Color Indicator
5. Test Strip Slot


5. To Set Up

1. Download and install the Pearl Fertility app on your mobile device. For a list of compatible devices, scan the QR Code at right:
2. Open the mobile app and register your device using the activation QR Code located in the kit's instruction sheet (included separately). Continue by following the on-screen instructions in the mobile app.
3. Note that you can only interpret test results using the mobile app.

6. How It Works

This system uses hormonal profiles, acquired by testing hormones FSH and LH (follicle stimulating hormone and luteinizing hormone) for a few days, to determine when ovulation is most likely. The system then displays the Fertile Window (the daily chance to conceive), based on those hormonal readings.

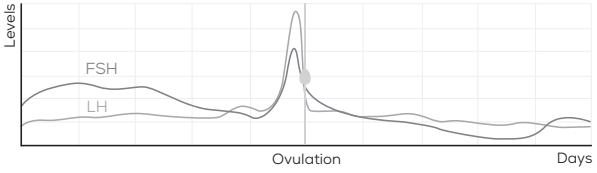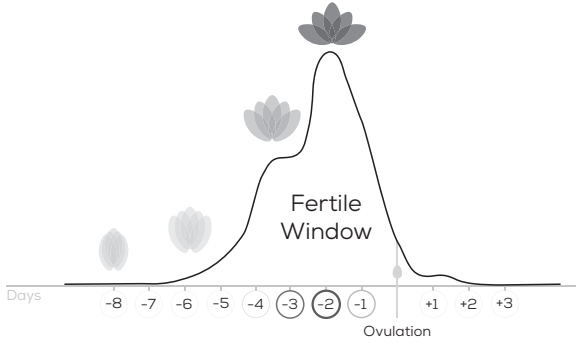[1] Treloar, A. E., et al (1967). Variation of the human menstrual cycle through reproductive life. Int J Fertil, 12(1 Pt 2), 77-126
[2] Stanford, et al "Timing intercourse to achieve pregnancy: current evidence." Obstetrics & Gyn. 100.6 (2002): 1333-134

7. Symbols Used

	Ovulation
	Intercourse
	Protected Intercourse
	More information button
	LH test strip (pink)
	Test completed (pink)
	Today button

	Heavy flow/ Cycle start
	Medium flow
	Light flow
	Spotting
	FSH test strip (blue)
	Test completed (blue)
	Pregnancy test strip (green)

	Expose the test strip to the first morning urine

	time duration

Questions?
The Pearl team is on hand to help! For more information, visit www.pearl-fertility.com/faq. The mobile app also has a contact function in your profile to send us questions and feedback directly. Our website and the information sections in the app offer some answers and background information about your cycle.


8. To Use

The Pearl Fertility App has three main screens: Calendar view, Cycle view, and Measurements view. Calendar view displays cycle information and the Fertile Window in calendar format so you have all important information at a glance. Cycle view shows your current cycle, so you are in control of current events. Measurements view shows how the tests are interpreted, so you gain insights every day. Before using the app, it is important to know which test strips to use each morning. The test strips need to be exposed to first morning urine. Instructions for using the test strips are included in a separate instruction manual.

1. **Wait for your cycle to start.** On the first day of your cycle (the first day you notice bleeding), tap on the flower to notify the app that your period has started. To do this, select a drop on the Period Input menu; the app will ask you if this day marks the beginning of your cycle. Select yes. The Cycle view will refresh starting from Day 1.


2. **Start testing when you wake up the next morning.** Starting on Day 2 of your cycle (the day after you first notice bleeding), the app will tell you which hormone tests you need to use that day. The symbols on the Cycle View days indicate which test strip to use (pink for LH and blue for FSH). Always use the first urine of the morning and read the instructions for each test strip.

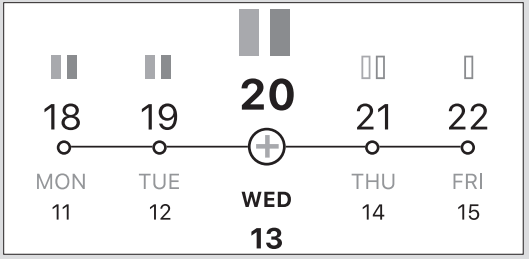
You have two different testing methods: using the strip holder and using a urine container.

To use the strip holder:

1. Slide the Protective Cover back completely and place the appropriate test strip(s) (LH and/or FSH) into the color-coded strip holder slot(s) all the way to the “Max” line. Make sure the “Max” marking aligns with the strip holder.

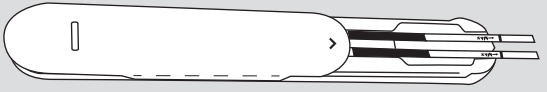2. Slide the Protective Cover back in place.

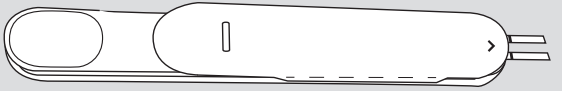3. Hold the tip of the test strip, facing downward, in the urine stream for about three seconds. Make sure the arrow marking on the strip holder does not get wet.

Clean the strip holder after each use (for more information see chapter on cleaning).

To use a urine container:

1. Collect your morning's urine sample in a clean, dry container.
2. Dip the tip of the test strip(s) into the urine up to the "Max" line for about five seconds.


After use:


1. Place the test strip(s) on a non-absorbent surface such as the blank side of the test strip packaging.
2. Use the timer function in the mobile app to measure the evaluation time. Wait 10-15 minutes for the test strips to develop (but no longer than 20 minutes). After 20 minutes, results are not accurate and the test will need to be redone.

Interpret the results:

After testing, select the appropriate strip(s) in the Input Menu on the mobile app. Take a photo of the test strip by centering it in the rectangle of the app window. Make sure the strip is well lit and is free of shadows and glare.

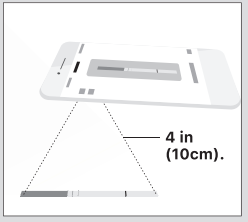or

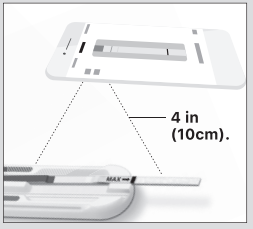Align the colored section of the test strip with the colored section in the app or let the app's Automatic Strip Detection feature find it for you. You can pinch, drag, and rotate the image to ensure it is captured correctly.


Drag the control handles on the image to align the test lines with the control lines in the mobile app. Try to center them as much as possible.


After aligning the lines, click “Done” to finish; the strip image will be saved and a tick mark will display on the input dialog, indicating you have completed the day’s test.

Dispose the test strip(s) by sliding the Protective Cover back while holding the Strip Holder vertically. The test strip(s) should come out of the holder; tap on it if the test strip does not come out easily.

As you continue testing, the mobile app will develop your hormonal profile and update your Cycle View in the app. The flower symbol represents your fertility chance for that day.

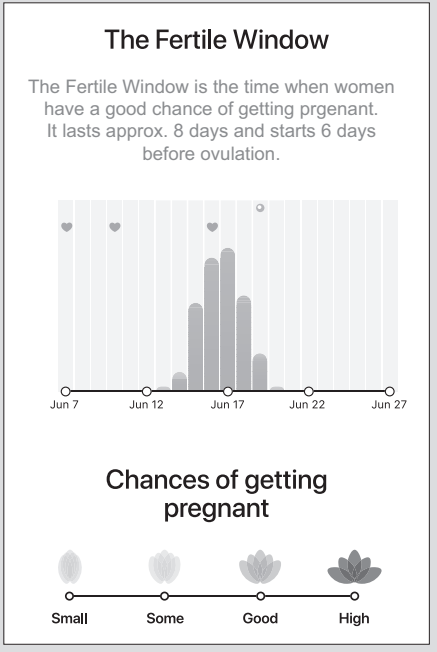
The term “fertile days” includes the days within the cycle when there is a good chance of pregnancy.This periode lasts approximately 8 days - it begins about 6 days before ovulation and ends one day after ovulation. The kit can tell you up to 6 days in advance when is your ovulation, covering up to 8 days of your fertile window (see chart).
The opening of the petals represents the probability of getting pregnant as shown in the Fertile Window. The calendar colors will help you understand how this Fertile Window works.

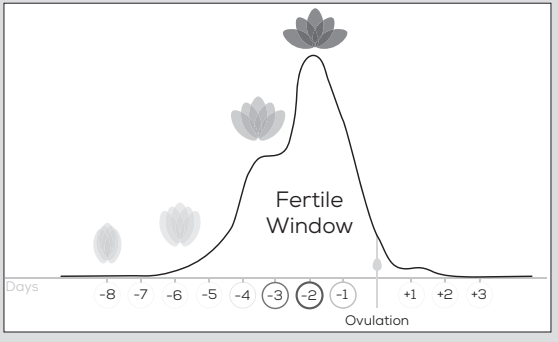9. Care, Maintenance and Disposal

1. Slide the Protective Cover by pulling it and the Strip Holder in opposite directions.

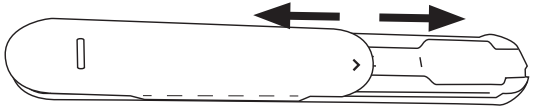2. Clean the strip holder with water after each use.
3. Dry the strip holder with a soft, dry cloth. Do not disinfect the Strip Holder or use any cleaners, solvents, or boiling water.
4. Slide the protective cover (with the arrow mark facing forward) back onto the Strip Holder.

Dispose of the test strips and the Strip Holder in your normal household waste after your cycle is complete.

Disposal

Observe the local regulations for material disposal. Dispose of the device in accordance with local regulations.
If you have any questions, please contact the local authorities responsible for waste disposal.

10. Technical Specifications

Model	OT 80
Type	C-BEU
Storage and transportation conditions	Temperature 39.2°F - 131°F (4°C – 30°C)
System requirements iOS ≥ 11.0, Android™ ≥ 6.0, List of compatible devices:	

For more information regarding our product line in the USA, please visit: www.beurer.com

Distributed by: Beurer North America LP
1 Oakwood Boulevard, Suite 255, Hollywood, FL 33020, USA
www.beurer.com

Questions or comments? Call our US-based customer service toll free at 1-800-536-0366.

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