


BELLAVITA BATHLIFT NOTES ON THE CORRECT USE OF THE BACKREST

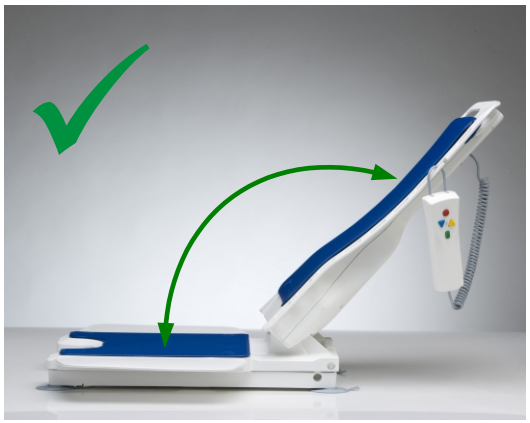


Ensure the bathlift base is fully lowered before AND during any operation on the backrest ...

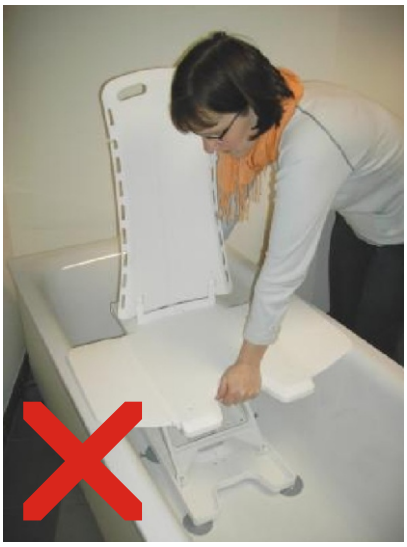

Do not install, remove or alter the backrest unless the bathlift is fully lowered and reclined ....

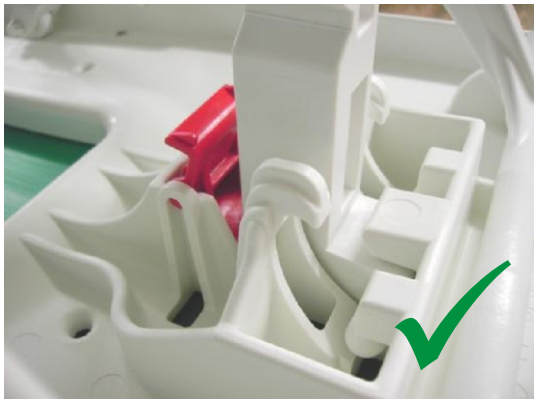

... as this is to ensure the gear rack locates properly in the ramp section (shown above) ...

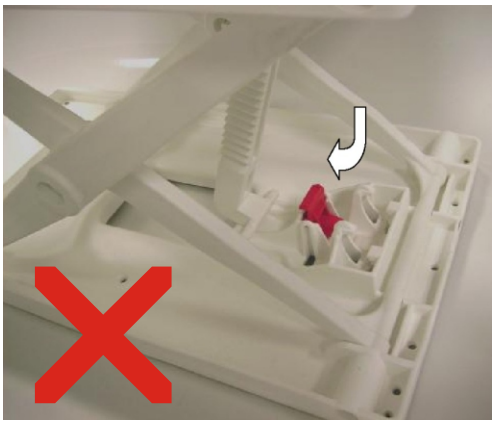

... or this may cause the gear rack to incorrectly locate ...



... and the bathlift will raise with an almost upright back and no risk of damage



... and the backrest to raise at an angle, risking damage to the device.