
*Single Sided Printing

Front of the QSG


50mm

75mm

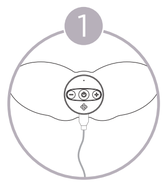


1

Unpack device, connect type-C cable and charge fully before use.

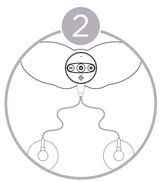


2

Attach cord with two electrode pads and clip unit to butterfly pad.

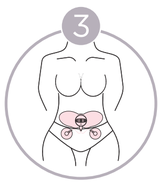


3

Remove protective film and attach gel pads to areas of treatment.

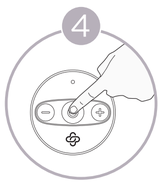


4

Hold down the power button for 2 seconds to turn it on.

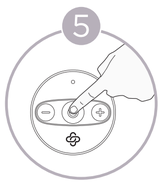


5

Press the power button to choose your desired treatment mode.

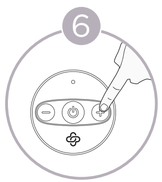

6

Use the + or - button to adjust the intensity and start the treatment.

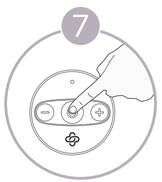

7

Press and hold the power button for 2 seconds to switch off.

Scan for support

WW_WW-300_UK_QSG_20250501_v2

54mm

79mm

kinetik WELLBEING

COMPLETE PERIOD PAIN RELIEVER

WITH TENS TECHNOLOGY

SKU WW-300

REF AD-2300

Cover

54mm

79mm

kinetik WELLBEING

For more information visit:
kinetikwellbeing.com

Back Cover