Solo Bed Transfer Aid



VY428/VY438/VY428N/VY438N

Fixing and Maintenance Instructions

EN

Downloaded from www.Manualslib.com manuals search engine

INTRODUCTION

Thank you for deciding to purchase the Aidapt Solo Bed Transfer Aid. The Bed Transfer Aid is manufactured from the finest material and components available, when installed and used correctly it is designed to give many years of reliable trouble-free service.

The Solo Bed Transfer Aid can be fitted to single, double, queen and king size beds. It is available in two finishes; white epoxy and chrome.

The Solo Bed Transfer Aid is not designed to prevent users falling from their bed.

PLEASE NOTE: VY428/VY428N IS THE SOLO BED TRANSFER AID WITHOUT STRAP

NB. This equipment must be installed by a competent person and a risk assessment undertaken as to the suitability of the Bed Transfer Aid for the particular user and user's bed. We would refer customers to 'Bed Rail Risk Management' LAC Number 79/8 published by the HSE.

VY428/VY438	**127 kg (20 st.) weight limit**
VY428N/VY438N	**159 kg (25 st.) weight limit**


Note: Do not exceed the weight limit stated – doing so could put the user at risk.

BEFORE USE

- Carefully remove all packaging. Avoid using any knives or other sharp instruments as this may damage the surface of the product.
- Check the product for any visible damage If you see any damage or suspect a fault, do not use your product as it may compromise your safety, but contact your supplier for support.

ASSEMBLY AND INSTRUCTIONS

1. Your Bed Transfer Aid is shipped ready assembled with its handle at its lowest setting (VY428 & VY438).
2. To install, remove the mattress from the bed where the Bed Transfer Aid is to be installed and check for any obstructions.
3. Place the Bed Transfer Aid onto the base of the bed with the large flat base on the bed frame and the looped handle facing upwards level with the side of the bed.
4. Take the shorter of the two straps and attach to the shorter horizontal crossbar level with the side of the bed by feeding the female buckle through the large plastic loop and pulling tight (see fig. 1).

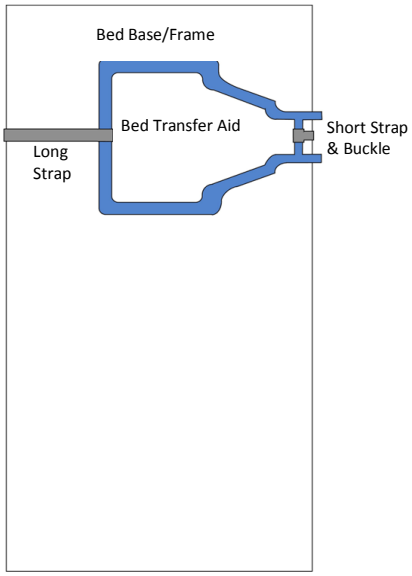


Downloaded from www.Manualslib.com manuals search engine

5. Take the longer of the two straps, remove the male buckle taking care to note the correct method for reattachment, then loop the longer strap around the longer horizontal crossbar situated in the centre of the bed frame and pull tight (see fig. 2).
6. Reattach the male buckle to the longer of the straps and feed the buckle and strap underneath the base of the bed (see fig. 3).
7. Move the Bed Transfer Aid up and down the bed to find the best position for the user to safely lift themselves.
8. Ensure that the frame of the Bed Transfer Aid is square to the side of the bed and that both of the black feet at the base of the unit are touching the side of the bed.
9. Clip the two buckles together and tighten the strap ensuring the Bed Transfer Aid is secure in its position and the buckle is clipped securely. Tuck away any excess strap to avoid trips.
10. Replace the mattress ensuring that the Bed Transfer Aid handle is pointing upwards and the unit is in the optimum position to suit the client's needs.
11. If required a second strap set can be fitted (VS247).


HEIGHT ADJUSTMENT (VY428 & VY438)

Height of the handle on the Bed Transfer Aid can be adjusted to suit the client's requirements. To adjust the handle simply remove the two grey e-clips from the base of the handle and lift the required height ensuring that you do not exceed the last pre-punched hole on the base of the frame.

Once the desired height is achieved, re-insert the e-clips into the pre-punched holes, ensuring that both clips pass completely through the tube frame, the inner tube and out the other side and are in place securely as shown below.

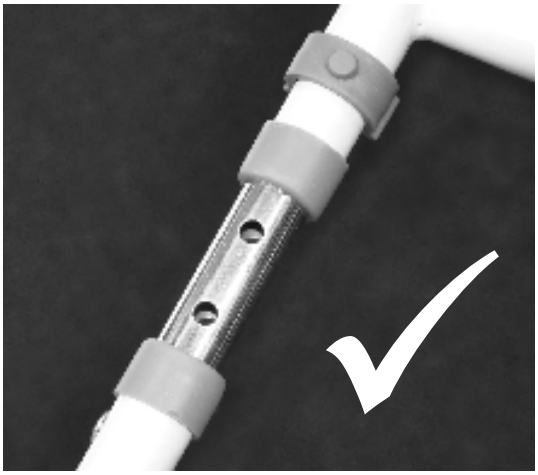

CORRECT

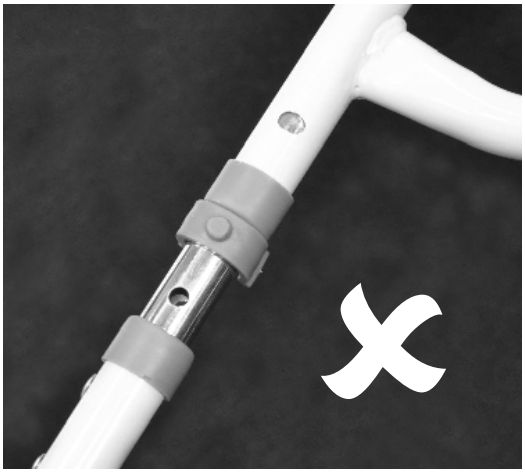

INCORRECT

Check that the unit is properly adjusted by pulling the handles upwards, if there is no movement, the height adjustment has been carried out correctly.

Downloaded from www.Manualslib.com manuals search engine



Fig. 1



Fig. 2



Fig. 3

CLEANING

The Solo Bed Transfer Aid should be cleaned using a non-abrasive cleaner or mild detergent with a soft cloth. Only use a mild detergent or disinfectant; never use an abrasive cleaner. Do not exceed 50ºC. Always ensure that you wipe dry all equipment and treat joints with a water dispersing agent, e.g. WD40.

REISSUE

If you reissue or are about to reissue this product, please thoroughly check all components for their safety.

This includes:

- Tightness of all nuts/bolts/castors and other screw in/bolt in/push in components
- Upholstery (if applicable) for security, splits, exposed staples, etc
- Straightness and integrity of all tubing and welds

If in any doubt, please do not issue or use, but immediately contact your supplier for service support.

When reissuing, a risk assessment must be undertaken as to the suitability of the Bed Transfer Aid for the particular user and user's bed.

IMPORTANT INFORMATION

The Bed Transfer Aid hand grips have been designed to meet the latest MHRA guidelines on entrapment, however you must ensure that there is no risk to the user or carer after the installation and service of the Bed Transfer Aid with Strap or without Strap.

Do not allow children or unauthorised persons to play with the equipment or use the bed transfer aid without proper supervision.

CARE & MAINTENANCE

Please make a safety check of the product at regular intervals or if you have any concerns.

Downloaded from www.Manualslib.com manuals search engine