Fixing and Maintenance Instructions
Shower Chairs

N.B. This equipment must be installed by a competent person.

Instructions for:

VB499/VB500	Milton Shower Chair	190 kg (30 st.) weight limit
VB501	Luton Shower Chair	190 kg (30 st.) weight limit
VB648/VB506/VB650	Hythe Shower Chair	190 kg (30 st.) weight limit
VB507	Charing Shower Chair	190 kg (30 st.) weight limit
VB510	Darenth Shower Chair	190 kg (30 st.) weight limit
VB511	Round Adjustable Shower Stool	115 kg (18 st) weight limit
VB509	Rochester Shower Stool	190 kg (30 st.) weight limit
VB508	Dartford Shower Stool	190 kg (30 st.) weight limit
VB635	Strood Shower Stool	159 kg (25 st.) weight limit
VB636	Pembury Corner Shower Chair	159 kg (25 st.) weight limit

Do not exceed the weight limit stated – doing so could put the user at risk.

1. All Shower Chairs/ Stools

Please examine your shower chair for any visible damage before using. If you see any damage or suspect a fault, please contact your supplier.

In this case, please do not use your product as it may compromise your safety.

2. Cleaning

All of the product components are manufactured from quality materials. However, it is very important that you only clean your shower chairs and stools using a non-abrasive cleaner or mild detergent with a soft cloth avoid using lime scale remover. Abrasive cleaners e.g. AJAX and or abrasive cleaning pads will seriously damage your product beyond repair.

Certain products have clip on/clip off backs and/or seats, you can easily remove these if required, clean with a mild detergent and soft cloth in a warm water solution. Do not exceed 50°C when cleaning the seat cover and/or back as this will permanently damage the plastic.

N.B. Regular daily cleaning will reduce the build up of limescale on your product.

Under No Circumstances use any abrasive cleaners or materials on any of the components. Also check that the cleaner you intend to use is suitable for ABS, Vinyl, Polymer Plastics and Thermo plastic rubbers and Polyurethane Rubber.

3. Adjustments to Extension Legs

Remove the 'E' clip from the extension leg and extend all four legs evenly to the desired height ensuring you do not exceed the last punched hole of the extension leg. (see fig. 1)

Once you have achieved the desired height; reinstall the 'E' clip, ensuring that it has passed through the extension leg and is fitting neatly and snug around the outside leg. (see fig. 2).

Never install the extension legs unevenly or slanting, this will compromise the safety of the user.

Aidapt Bathrooms Ltd., Lancots Lane, Sutton, St. Helens, Merseyside, WA9 3EX
Telephone: +44 (0) 1744 745 020 • Fax: +44 (0) 1744 745 001 • Web: www.aidapt.co.uk
E-mail: accounts@aidapt.co.uk • adaptations@aidapt.co.uk • sales@aidapt.co.uk

4. Use of Bath Oils

Caution should be observed if you are using any form of bath oils or oily skin preparation with your shower chair as this may impair your natural contact skin resistance.

In the event of using any such oil, your seat must be thoroughly cleaned with a suitable detergent to break down any oily residue. Failure to observe this recommendation could impair your safety.

5. Reissue

If you reissue or are about to reissue this product, please thoroughly check all components for their safety.

This includes:

- PAT test for electrical items
- Tightness of all nuts/ bolts/castors
- Other screw in/ bolt in/ push in components.
- Also check all upholstery for security, splits, etc.

<u>If in any doubt, please do not issue or use, but immediately contact your supplier for service support.</u>

fig.1


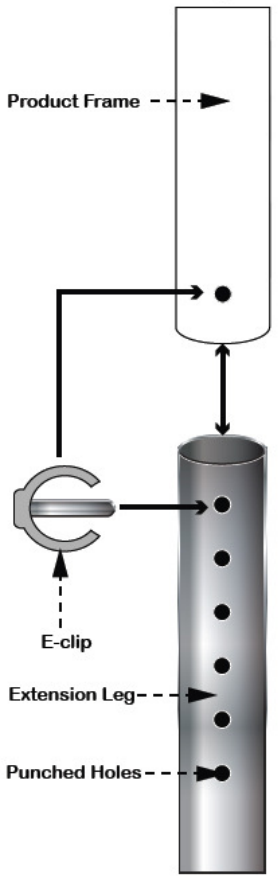


fig.2

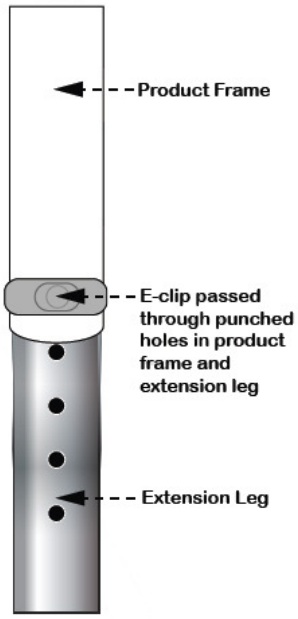

6. Maintenance

Please make a safety check of the product at regular intervals or if you have any concerns.

Aidapt Bathrooms Ltd., Lancots Lane, Sutton, St. Helens, Merseyside, WA9 3EX
Telephone: +44 (0) 1744 745 020 • Fax: +44 (0) 1744 745 001 • Web: www.aidapt.co.uk
E-mail: accounts@aidapt.co.uk • adaptations@aidapt.co.uk • sales@aidapt.co.uk